



3100 WristOx[®]


Wearable Digital
Pulse Oximeter

Wearable Oximetry for the Ambulatory Patient

NONIN's WristOx[®] 3100 is a small, lightweight pulse oximeter designed to be worn comfortably on the patient's wrist. Ideal for monitoring daily activities and for overnight studies, the compact design of the WristOx is intended to give the patient greater freedom while improving patient compliance. With impressive battery life, powerful memory, and NONIN's proven accuracy - the WristOx provides a cost effective solution for both long and short term monitoring. Simply slip the WristOx onto your patient's wrist and begin monitoring SpO₂ and Pulse Rate instantly. Right out of the box in the spot-check mode, the WristOx is a powerful tool: to complete the picture, add nVISION[®] software and enable advanced WristOx features.

Product Highlights

- ✓ Small - Fits comfortably on the patient's wrist
- ✓ Versatile - Ideal for short or long term monitoring applications
- ✓ Powerful - 24 hours of battery life and 33 hours of non-volatile memory
- ✓ Easy-to-Use - Simply plug in the sensor and go; large easy-to-read LCD
- ✓ Expandable - Advanced features available when used with nVISION
- ✓ Durable - 3-year warranty
- ✓ Compatible - Works with NONIN's full line of PureLight[®] sensors

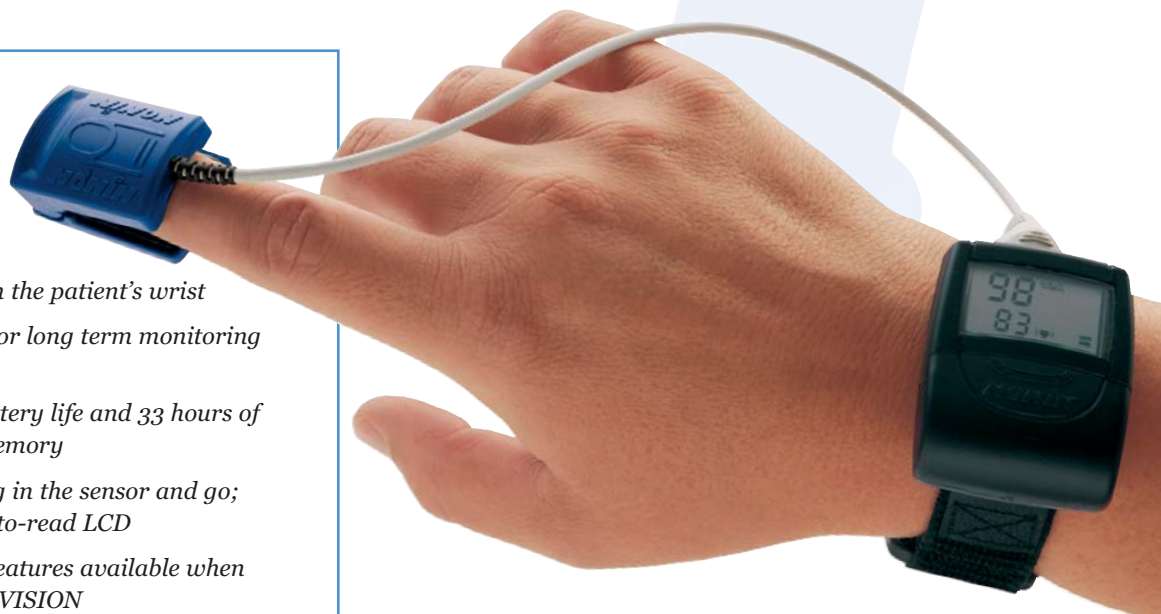

Versatility

The WristOx® 3100 is an ideal choice for multiple applications.

Hospital

Cardiopulmonary Rehabilitation - Six-Minute Walk Assessment,
Long Term Oxygen Therapy (LTOT) Monitoring
Sleep Labs - Pre and post PSG Studies, Adult & Pediatric pre-surgical assessments
Step Down Units - Ambulatory Monitoring

Physician Office

LTOT Monitoring- Oxygen
Certification/Recertification,
Oxygen Titration, Overnight Studies,
Pre-Surgical Assessments

Nursing Home

LTOT Monitoring - Patient
Assessment, Oxygen
Certification/Recertification,
Oxygen Titration

Advanced Features with nVISION®



Combining the WristOx with NONIN's nVISION® Data Management Software will unleash the full power of your wrist oximeter. nVISION is a Windows®-based software application which provides data analysis and report generation when used with NONIN's memory products. Purchase nVISION with your WristOx and enjoy advanced features such as:

- ✓ Sensor Activation Mode (*sensor connect/disconnect turns WristOx On/Off*)
- ✓ Programmed Mode (*for pre-programming up to 3 separate Start/Stop times*)
- ✓ Set Date, Time and Patient ID
- ✓ Set Patient Data Storage Rate (*for choosing trend resolution*)
- ✓ Choose a Full or Partial Display (*turns off SpO₂ and Pulse Rate display for anxious patients*)
- ✓ Access 33 Hours of Stored Data (*for sleep studies or ambulatory monitoring*)
- ✓ Clear Memory

Select Your Options

- ✓ 8000J-WO Adult Flex Sensor
- ✓ N-cell Batteries (*box of 24*)
- ✓ Reusable Wristbands (*pack of 10*)
- ✓ 5-inch Reusable Extender Bands (*pack of 10*)
- ✓ Adapter Cable (*for compatibility with all NONIN PureLight® sensors*)
- ✓ nVISION Data Management Software
- ✓ PC Interface Cable
- ✓ WristOx Carrying Case



Please contact your authorized distributor for more information.

Specifications

Oximeter

Oxygen Saturation Range (%SpO ₂)	0% to 100%
Pulse Rate Range	18 to 300 pulses per minute
Displays	
Numeric Displays	3-digit Indicators
Pulse Indicator	Pulse Strength Bargraph
Accuracy	
Blood Oxygen Saturation (%SpO ₂) (±1 SD)*	70 - 100% ± 2 digits
Pulse Rate	±3%
Temperature	
Operating	+32° to +122°F (+0° to +50°C)
Storage/Transportation	+14° to +122°F (-10° to +50°C)
Humidity	
Operating	10% to 90% noncondensing
Storage/Transportation	10% to 95% noncondensing
Altitude	
Operating Altitude	Up to 40,000 feet
Hyperbaric Pressure	Up to 4 atmospheres
Power Requirements	Two 1.5V alkaline N-cell batteries
Battery Life	
Operating	Minimum 24 hours of continuous operation
Storage	9 months
Dimensions (without sensor or strap)	1.75" wide x 2" high x .75" deep 4.445 cm x 5.08 cm x 1.905 cm
Weight (without batteries or wrist strap)	~.88 ounces (24.95 g)
Memory	
Type	Nonvolatile
Capacity	33 hours @ 4 sec. data storage rate 16 hours @ 2 sec. data storage rate 8 hours @ 1 sec. data storage rate
Warranty	3 years

Classifications per IEC 60601-1/CSA601.1/UL2601-1

Type of Protection	Internally powered (on battery power)
Degree of Protection	Type BF-Applied Part
Mode of Operation	Continuous

* +/- 1 Standard Deviation encompasses 68% of the population.

Quality systems registered to ISO 13485:2003 and EC certificates per Annex II, Clause 3 of EC Directive No. 93/42/EEC concerning medical devices.



Specifications are subject to change without notice.



Your Partner of Choice



With over 20 years of experience and dedication to the design and support of noninvasive monitoring devices, NONIN has helped many medical professionals meet clinical and economic objectives. Trusted by clinicians worldwide, countless tests are performed every day on thousands of NONIN pulse oximeters in more than 125 countries. NONIN's dedication to technological leadership, precision manufacturing and uncompromised customer support ensures quality products and service you expect.



Nonin Medical, Inc.
13700 1st Avenue North
Plymouth, MN 55441-5443, USA
Phone: +1.763.553.9968 Toll Free: 1.800.356.8874
Fax: +1.763.553.7807
info@nonin.com