

# AED Plus® Trainer2

## Training Accessories



### Complete Simulation Control at the Instructor's Fingertips

The ZOLL® AED Plus Trainer2 allows CPR and first aid instructors to have students experience what it's like to save a life with the AED Plus. Choose from four different fixed rescue scenarios or have complete manual control of the entire rescue. A wireless remote lets you control multiple AED Plus Trainer2s.

With this simple remote, instructors can set up a rescue scenario and download to the AED Plus Trainer2 in seconds. The student pushes the on-off button on the Trainer2 unit to begin. The instructor then determines whether to shock on the next simulated heart analysis or not, whether the pads have been properly attached or not, and when the Trainer2 unit should say "PUSH HARDER" and "GOOD COMPRESSIONS." A simple pause button lets the instructor stop the rescue altogether, confer with the student, and then pick up exactly where the simulated rescue left off.

### 2010 AHA and ERC Guidelines Ready

Trainer2 software easily switches from AHA to ERC 2010 Guidelines

### The AED Plus Trainer2 Includes Everything an Instructor Needs

- AED Plus Trainer2 training AED
- Wireless remote controller
- One set of AED Plus Training CPR-D-padz®
- One set of AED Plus Trainer CPR-D-padz reusable gels
- Four D cell batteries
- Two AA cell batteries
- Operator's Guide

## AED Plus Trainer2 Accessories



AED Plus Trainer2 Replacement Remote  
Part# 8008-0007



AED Plus Trainer Carry Case  
Part# 8000-0375-01



AED Plus Trainer Demo Manikin  
Part# 8000-0835-01

## AED Plus Trainer2 Specifications

**Dimensions:** 5.25"H x 9.5"W x 11.5"D  
(13.3 cm H x 24.1 cm W x 29.2 cm D)


**Weight:** 4 lbs 14 oz (2.2 kg)

**Power:** 4 D cells (included) for Training AED  
(Good for more than 100 hours of training)  
2 AA cells (included) for remote controller  
AC adapter not available—battery operated only.



## Rescue Scenarios

Instructors can set up five different rescue scenarios with this simple remote controller.

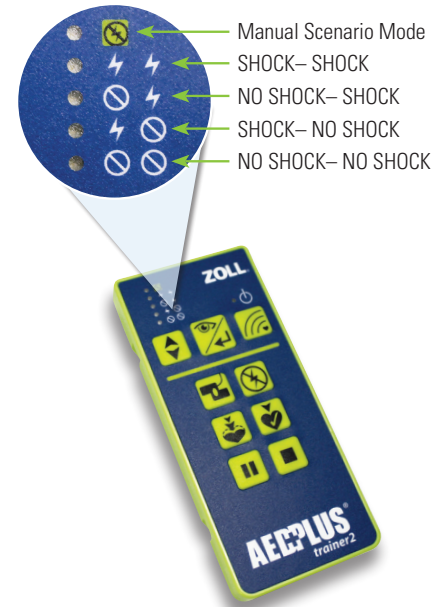

## AED Plus Trainer2 Accessories

(Replacements must be ordered separately)



Replacement Adhesive Gels for  
CPR-D-padz Trainer Electrodes  
Part# 8900-0803-01



CPR-D-padz Trainer Electrodes  
Part# 8900-0804-01



Stat-padz® II Trainer Electrodes  
Part# 8900-0805-01  
(case of 6)



Pedi-padz® II Trainer Electrodes  
Part# 8900-000861-01  
(case of 6)

**ADVANCING  
RESUSCITATION.  
TODAY.®**

ZOLL Medical Corporation  
Worldwide Headquarters  
269 Mill Road  
Chelmsford, MA 01824  
978-421-9655  
800-348-9011

For subsidiary addresses and fax numbers,  
as well as other global locations, please  
go to [www.zoll.com/contacts](http://www.zoll.com/contacts).



Specifications subject to change without notice.

©2011 ZOLL Medical Corporation. All rights reserved. Advancing Resuscitation Today, AED Plus, CPR-D-padz, Pedi-padz, Stat-padz, and ZOLL are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.

Printed in U.S.A. 061110 9656-0202

