



## ELITE-LINE MULTI-USE IMAGING POWER TABLE with 3-Section Top & Drop Window



ADI903-04



**Multi-Use Imaging Power Table  
with 3-Section Top & Drop Window**

Level up your imaging equipment with our Viva Comfort Multi-Use Imaging Power Table with a Three-Section Top & Drop Window – your premier solution for versatile and comfortable imaging procedures. The three-section top design includes a pneumatic adjustable backrest, manual footrest, and a convenient drop window, providing optimal patient positioning for various procedures. To give you a premium imaging table that truly is unmatched, we covered the outfitted 3" padding with antimicrobial upholstery. With an included easy-to-use hand control, the table's height adjustments are effortlessly accessible. The imaging table's design enables a swift and smooth transition in just 16 seconds for height adjustments without load and 27 seconds with a remarkable weight capacity of 550 lbs. The table's adaptable features extend to its Trendelenburg and Reverse Trendelenburg angles, as well as the leg section and backrest adjustments, guaranteeing patient comfort during various medical procedures.

We built the power table specifically with precision and stability in mind, and that is why our base and top frame are all-welded and equipped with a dual-column lift system that solidifies its exceptional rigidity. We outfitted the steel base with 3" heavy-duty wheels to withstand consistent use, handle the total weight capacity, and allow for swift maneuverability. As a bonus, you can buy add-on accessories that will customize the table according to your specific clinical requirements.

## FEATURES

**Design:** Three-section top design with pneumatic adjustable backrest, manual footrest, and a convenient drop window.

**Sanitation:** Antimicrobial upholstery that is fire, stain, mildew, oil, water, and wear-resistant.

**Adjustable:** Easy-to-use hand control for height adjustments ranging from 21.5" to 37.5".

**Adaptable:** Offers Trendelenburg and Reverse Trendelenburg angles of +20° to -20°, along with leg section and backrest adjustments of up to 80°.

**Stability:** All-welded base and top frame, equipped with a dual-column lift system, provide exceptional rigidity, precision, and solidity.

**Optional Add-On's:** Additional accessories like our side rails, stirrups, paper roll holder and cutter, and safety straps are available for purchase, allowing customization according to clinical requirements.

# SPECIFICATIONS

## Elite-Line Multi-Use Imaging Power Table with 3-Section Top & Drop Window

<b>Material</b>	Metal, Leather, High-Density Foam, ABS Caster
<b>Overall Dimensions</b>	21.45"-37.2" H x 27.56" W x 76.77" L
<b>Backrest Dimensions</b>	36.6" x 27.56"
<b>Footrest Dimensions</b>	15.75" x 18.9"
<b>Mid-Section Dimensions</b>	21.26" x 27.56"
<b>Window Dimensions</b>	10" x 10"
<b>Net Weight</b>	242.51 lb.
<b>Weight Capacity</b>	550 lb.
<b>Voltage</b>	110v
<b>Watts</b>	400VA
<b>Motor Certification</b>	CE, ETL
<b>Padding</b>	3" thick high-density foam
<b>Upholstery</b>	Antimicrobial
<b>Color Options</b>	Black and Midnight Express
<b>Power Cord Length</b>	13.12 feet

## FEATURES



Adjustable Backrest



Hand Control



Drop Window



Cable Management




Heavy-Duty Wheels



Dual-Column Lift

## OPTIONAL ADD-ON ACCESSORIES

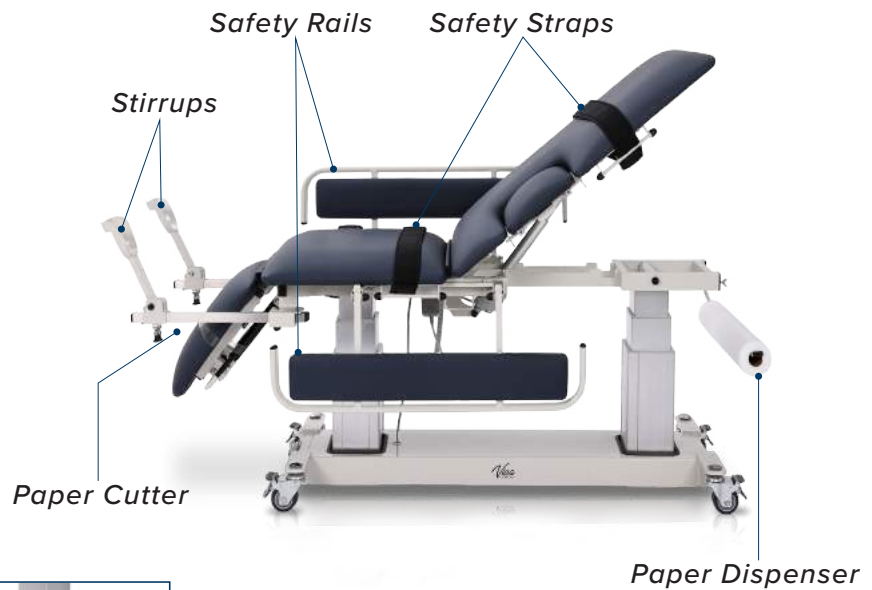

## UPHOLSTERY COLOR OPTIONS



Black



Midnight  
Express