

CRYODOSE™

by Nuance Medical



KEY FEATURES

- -70°C
- Non-flammable, ozone safe, no hazmat
- Extended shelf life
- Dosing cones with splash guards
- Replacement products available
- Pleasing look for greater patient acceptance



UNIQUE TO CRYODOSE

- 1 Tapered nozzle channel for secure fit
- 2 Responsive trigger to control discharge and avert waste
- 3 Six dosing cones with measurement markers
- 4 Arrow and round buds for more treatment options
- 5 Handy quick-reference usage guides
- 6 Reusable Treatment Kit for ease of use and lower cost

The CryoDose™ system makes your cryotherapy practice more perfect.

Nuance Medical, LLC
5931 Sea Lion Place, Suite 113
Carlsbad, CA 92010

760 585-4849
info@CryoDose.com

Your Local Rep Name Here
Business Name, If Needed
Street Address
Suite
City or Town, ST, zip code
phone num-ber
email@address

CryoDose™ offers 360° benefits


for every medical professional and their patients

EASY USAGE

USING TREATMENT BUDS


Step 1

Gently squeeze trigger to prime the bud.

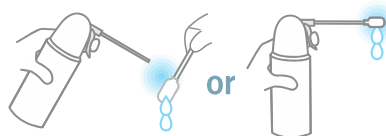

Step 2

Point bud downward 15-20 seconds to reach optimal freezing temperature.

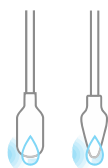

Step 3

Place bud onto lesion for recommended treatment time.


USING DOSAGE CONES

Step 1

Select appropriate cone to completely cover lesion.

Step 2

Fill between bottom 1/8"-1/4" fill lines.

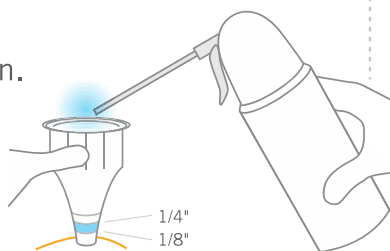

Step 3

Keep in place until evaporation is complete.

EASY ORDERING OPTIONS

CRYODOSE REUSABLE TREATMENT KITS

1101 162mL Kit w/ 6 cones, 40 buds

1103 162mL Kit w/ 80 buds

1501 236mL Kit w/ 6 cones, 40 buds

CRYODOSE CANISTERS

1100 162mL canister

1500 236mL canister

CRYODOSE BUDS

2001 mixed, 20 arrow, 20 round

2002 arrows, bag of 30

2003 round, bag of 30

CONES

3001 cones – 3, 5, 7, 9, 12, 16mm

*cryotherapy,
precisely*

EXCELLENT CPT REIMBURSEMENT

11200 \$85.13
skin tags

17000 \$81.25
benign lesions:
actinic keratosis, SK

17110 \$110.96
warts, plantar warts,
molluscum

54056 \$145.26
genital, penis, simple

56501 \$137.72
genital, vulva, simple

11201 \$19.40
each additional 10

17003 \$7.15
2 to 14

17111 \$131.39
15 or more

54065 \$232.64
penis, extensive

56515 \$227.32
vulva, extensive


CRYODOSE™

by Nuance Medical