




Closed Door System Fully Automatic Autoclaves

EZ9TM and **EZ10**TM

Autoclaves will meet all your sterilization needs

- Chamber door remains closed during entire cycle
- Hepa filtered air is circulated in chamber during drying
- All program parameters can be customized and stored
- New brighter and easier to read panel
- Two year parts and labor warranty

* Above photo shown with optional printer

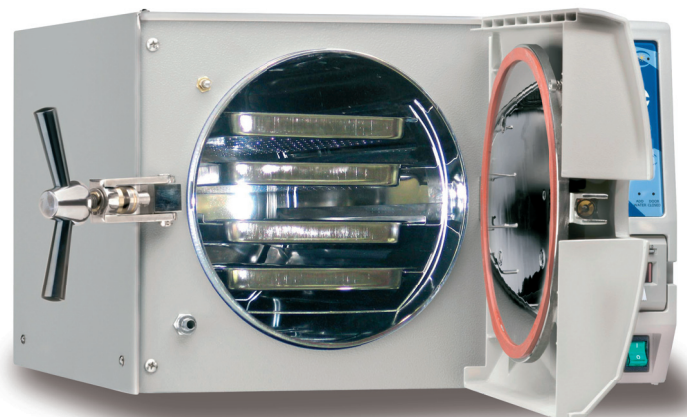

EZ10 with 4 large trays

With the simplicity of one touch design all your sterilization and drying needs are fulfilled. The Tuttnauer fully automatic autoclaves fill, sterilize, exhaust and dry at the touch of a button.

EZ9™ and EZ10™ Series:

| Item No. | Chamber Dimensions Ø x L In(mm) | Volume | Overall Dimensions DxWxH In(mm) | Standard Cassettes | Tray Dimensions DxWxH In(mm) | No. of Trays | Standard Unwrapped Cycle time | Voltage(V) Freq.(Hz) | Power(W) Current(A) | Shipping Weight lbs(kg) |
|----------------------|---------------------------------|----------|----------------------------------|--------------------|---------------------------------|--------------|-------------------------------|----------------------|---------------------|-------------------------|
| EZ9™ | 9" x 18" (230 x 470) | 5gal/19L | 21.5"x20"x14.4" (545x510x365) | 2 full/2 half | 16.3"x6.7"x0.8" (415x170x20) | 3 | Cold: 27 mins Hot: 13 mins | 120V 50/60 | 1400W 12A | 95 (43) |
| EZ10™ | 10" x 19" (254 x 475) | 6gal/23L | 21.5"x20"x14.4" (545x510x365) | 3 full/3 half | 16.3"x6.7"x0.8" (415x170x20) | 4 | Cold: 30 mins Hot: 14 mins | 120V 50/60 | 1400W 12A | 95 (43) |
| *EZ10k™ Kwiklave™ | 10" x 19" (254 x 475) | 6gal/23L | 21.5"x20"x14.4" (545x510x365) | 3 full/3 half | 16.3"x6.7"x0.8" (415x170x20) | 4 | Cold: 21 mins Hot: 11 mins | 230V 50/60 | 2200W 10A | 95 (43) |

*Model EZ10k™ requires a voltage between 220V and 235V for proper operation. When applicable a buck/boost transformer (0.5KVA) may be required.

Features:

- Closed door active drying system to maintain sterility and ensure efficient drying of packs and pouches
- High efficiency air pump circulates air in chamber during drying cycle
- 0.2µm Hepa air filter (replaceable)
- New brighter and easier to read panel
- Program parameters can be customized and stored

Control Panel:



1. Program Selection
2. Detailed Digital Display
3. Programmable Keys
4. System Indicators


Optional Internal Printer



Documents:
date, time,
temperature,
pressure.



Clean&Simple™ Ultrasonic Cleaners and Tablets

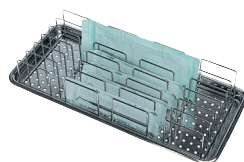

Instrument Pouch Rack



Steam Distiller



Chamber Brite Autoclave Cleaner