

Open Base Power XL Table with Adjustable Backrest

MODEL 80299

- Features easy-clean, 30" wide top with adjustable pneumatic backrest which is adjustable from either side and raises 80 degrees from flat
- 600 lbs. (272 kg) load capacity under normal use
- Smooth, electric, raise and lower height adjustment
- Pneumatic backrest cylinder is mounted internally
- Full width backrest hinge adds rigidity for greater support
- Heavy duty, steel frame is fabricated in the USA by Clinton Industries
- Dependable scissors-lift features an H-leg design for extra stability
- Rounded corner 3/4" plywood top with 2" of ultra firm padding covered by heavy weight upholstery
- Harden steel bolts and self lubricating pivot bearings for years of dependable service

CLINTON OPEN BASE & SHROUDED
POWERXL
Extra width, Extra capacity, Extra Affordable



80299 SPECIFICATIONS

Overall 72" L x 30" W x 18-35" H

Backrest 38" L x 30" W

- Low height may help achieve ADA tax credit

California residents please be advised of the following, pursuant to Proposition 65:

WARNING: This product can expose you to chemicals, including chromium, which is known to the State of California to cause cancer or reproductive toxicity. For more information go to www.P65Warnings.ca.gov.

1



Upholstery Colors



Slate Blue



Wedgewood



Royal Blue



Desert Tan



Burgundy



Country Mist



Warm Gray



Gunmetal



Black

OPTIONS

Positioning Straps Set

- 2" x 92" strap is secured with quick hook and loop Velcro® fasteners
- 2 straps per set (92" x 2")

080



Safety Rails (pair)

- 3 height positions
- Plunger height adjustment

082



Hand-Held Controller

- Alternative to foot control
- Stores easily under top when not in use

084



030 Paper Dispenser


- Accommodates paper rolls from 18" (45.72cm) to 21" (53.34cm) wide

030

Paper Cutter

- Clear plastic with elastic strap and buckle

040

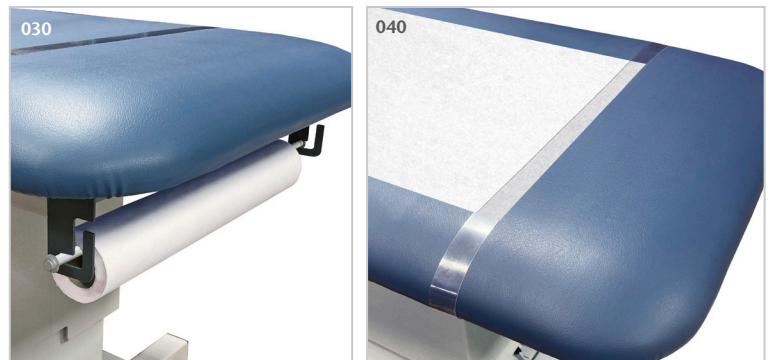

Paper Dispenser/Cutter Combo

- Paper dispenser and paper holder ordered as one option

034