



Tilting Power Hi-Lo Rehab Therapy Table



ADI903-05



The Tilting Power Hi-Lo Rehab Therapy Table is thoughtfully engineered to empower healthcare professionals in delivering precise, effective therapy to help diagnose conditions such as syncope, dysautonomia and physical and internal rehabilitation. This table is an essential tool for safeguarding patient health by facilitating treatments for conditions affecting the neurological and nervous system, heart, bladder, intestines, sweat glands, pupils, and blood vessels. Its electric Hi-Lo and tilt operation, controlled via a convenient handheld remote, offers seamless height adjustment from 22.4"–41" and tilt angles ranging from -2° to +90°. To guarantee we deliver optimal comfort to each patient, our thoughtfully crafted table boasts 2 inches of high-density foam, combined with an adjustable headrest section that adjusts from 0° to 75° supported by a reliable gas strut. Located at the end of the table is an outfitted footrest that can be easily adjusted, too, up to 20°. Each table includes an activity tray and removable chest, waist, and leg safety straps.

The therapy table supports up to 485 lbs. with the guide of two durable TiMOTION motors engineered to accommodate a broad spectrum of patient weights, providing reliability for diverse needs. Our Tilting Power Hi-Lo Rehab Therapy Table redefines reliability and excellence, making it an indispensable asset for every clinical setting to invest in.

Features

Outfitted: The integrated quick-release handle allows for the immediate flattening of the backrest to a supine position, providing a critical fail-safe for patients experiencing a sudden drop in blood supply or medical emergency.

Versatile: The table offers a height adjustment range of 22.4"–41" and tilt angles from -2° to +90°, ensuring it can adapt to various therapeutic applications.

Precise: An integrated inclinometer allows for accurate positioning, ensuring the table can be adjusted precisely to meet patient and doctor's needs.

Accessible: The 6.5" under-frame clearance ensures total compatibility with Hoyer lifts and other assistive devices, streamlining patient transfers and enhancing ease of use for the caregiver.

Comfortable: The table features 2 inches of high-density foam and an adjustable headrest that tilts from 0° to 75°, providing superior comfort for patients during extended use.

Functional: A removable activity tray is included, making it ideal for therapy exercises, fine motor skill development, and coordination or cognitive tasks.

Convenient: The table includes removable chest, waist, and leg safety straps, which can be neatly rolled and clipped to the sides when not in use for easy organization.

Durable: The 2-inch upholstery is crafted from antimicrobial PVC vinyl, which resists stains, mildew, oil, and wear, ensuring longevity and easy maintenance.

Safe: The design incorporates sleek, rounded frame legs to minimize injury risks and 5-inch lockable casters for effortless mobility and secure stabilization during use.

Robust: Powered by two TiMOTION motors, the table supports a weight capacity of 485 lbs., making it reliable for a wide range of patients.

Hygienic: A removable face cushion on the headrest section accommodates various treatments.

Efficient: The handheld remote makes adjustments easily, offering smooth transitions between height and tilt positions during therapy sessions.


Dimensions



Features



Handheld Remote

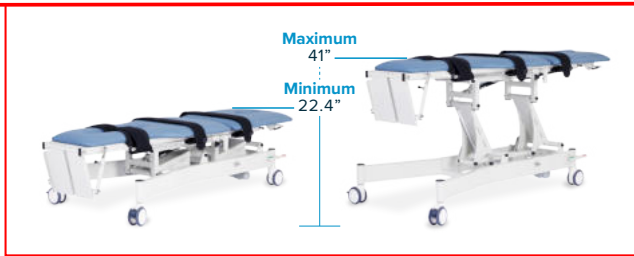

Adjustable Height



Adjustable Tilt Angle



Adjustable Headrest



Inclinometer



Removable Cushion



Footrest



Chest, Waist and Leg Security Straps



Activity Tray



5" Lockable Caster Wheels



Quick Release

Specifications

Materials	Steel Frame, Sponge, Leather, ABS Caster
Overall Dimensions	76.4" L x 25.2" W x 22.4-41"
Weight Capacity	485 lbs.
Net Weight	243 lbs.
Mattress Dimensions	76" L x 25" W
Head Section Dimensions	26" L x 25" W
Leg Section Dimensions	48" L x 25" W
Motors	2 TiMOTION Motors
Height Range	22.4"- 41" (By Handheld Remote)
Tilt Angle	-2° to +90° (By Handheld Remote)
Adjustable Backrest Angle	0 to 75° (By Gas Strut)
Footrest Tilt	0° to 20°
Footrest Size	12.5" W x 8" D
Wheels	5" Lockable Castors
Certification	CE, FDA, ISO 13485:2016
Included	Inclinometer, Removable Face Hole Pad, Activity Tray, Chest, Waist, and Leg Security Straps, Handheld Remote

Upholstery	Antimicrobial PVC Vinyl (Stain, Mildew, Oil, and Wear Resistant)
Color	Midnight Express
Padding Thickness	2"
Foam Density	36kg/m ³ (79 lbs/m ³)
Length of Cord	13'
Watts	380 VA
Voltage	100-240V
Decibels	Less than 65 dB
Activity Tray	29" L x 12" W
Waist Strap Dimensions	39" L x 4" W
Waist Strap Weight	0.65 lbs.
Chest Strap Dimensions	39" L x 4" W
Chest Strap Weight	0.65 lbs.
Foot Strap Dimensions	44" L x 4" W (The Widest Part: 8")
Foot Strap Weight	0.88 lbs.
Clearance From Floor To Table Base	6.5"