

PleurX catheter system

Vacuum bottle

Active vacuum technology drains quickly and comfortably without the need for gravity. Bottles are available in 500 mL and 1,000 mL sizes.



Roller clamp

The roller clamp provides easy control of drainage flow for optimal comfort.


Proprietary safety valve

The valve helps prevent inadvertent passage of air or fluid through the catheter. It also helps to reduce the risk of misconnection that may occur with a standard luer connection.



Polyester cuff

This promotes tissue ingrowth to help reduce infection risk and hold the catheter securely in place.

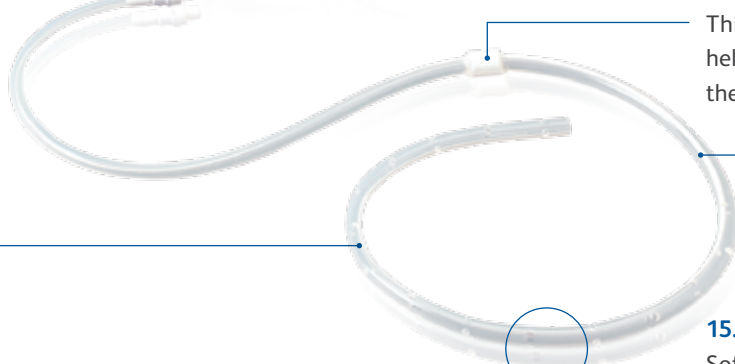


15.5 Fr Silicone catheter

Soft and flexible catheter conforms to the pleural or peritoneal space and minimizes insertion-site discomfort.

Fenestrations

Large, smooth fenestrations promote drainage and help prevent occlusions.

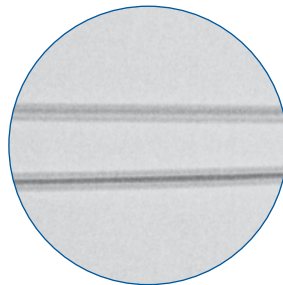

Barium stripe

Darker barium stripe provides more visibility under radiography.