

# BLICKMAN



## *Mayo Stands*

*Maximum utility in minimal space*

**Includes:**

- Chrome Mayo stands
- Stainless steel Mayo stands



# OUR MISSION

## **Customer Focus**

Exceed customer expectations through a commitment to continuous improvement of our products and services.

## **Quality Products**

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

## **Organizational Excellence**

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

## **Marketplace Recognition**

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

## **Collaborative Workplace**

Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.



**BLICKMAN**

# THE BLICKMAN DIFFERENCE

- Our Mayo stands have evolved to meet the needs of the marketplace for maximum utility in an increasingly crowded surgical environment.
- We have an array of options with varying material, bases, shapes, and adjustment mechanisms to suit your needs.
- All of our Mayo stands are mobile and use 2" roller blade, hard rubber, or twin wheel plastic hooded casters.

## Key Product Features

---

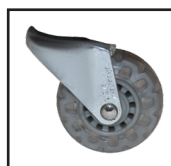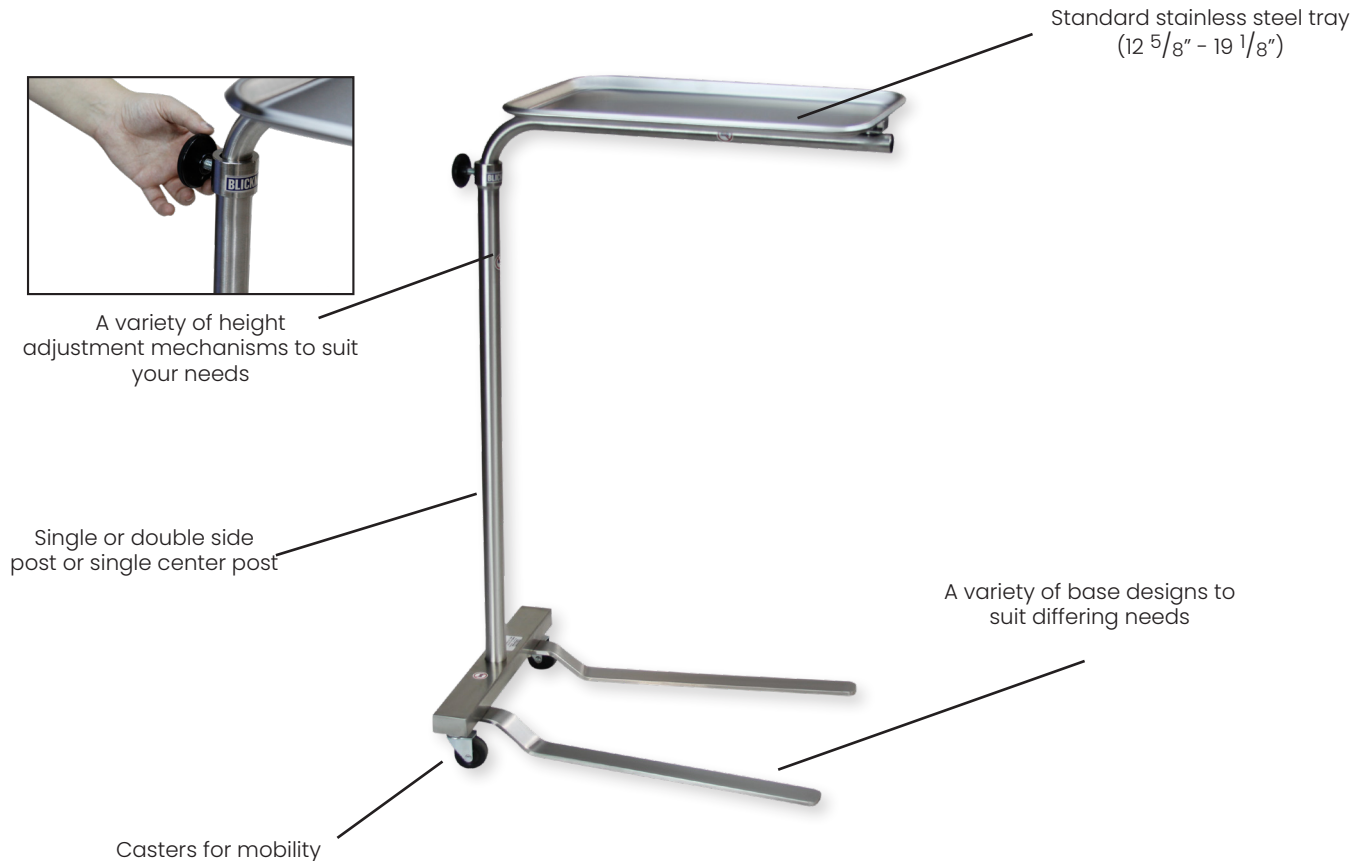

2" Rollerblade Casters



2" Hard Rubber Casters



2" Twin Wheel Hooded Casters

# Mayo Stands

Chrome

## CHROME MAYO STANDS

**Description:** Each of our chrome Mayo stands contain a 12 5/8" x 19 1/8" stainless steel tray and are height adjustable using the Tru-Loc™ friction knob. We offer different base and pole set ups to meet your varying needs.

### Middle Pole

**Base:** Five leg cast aluminum base with wall savers™ and 2" twin wheel hooded casters



### Double Pole

**Base:** Chrome loop base with two 2" twin wheel hooded back casters



### Single Pole

**Base:** Two pronged chrome base with single hard rubber back casters



Name	Middle Pole	Double Pole	Single Pole
Height Adjustment	29 1/2" - 49"	34 1/2" - 54"	34 1/2" - 54"
Model #	1501	1510	1515
Item #	0661501000	0661510000	0661515000
Shipping Weight	14 lbs.	18 lbs.	12 lbs.

# Mayo Stands Stainless Steel

## NEWARK MAYO STANDS

**Description:** Each of our stainless steel Mayo stands contain a 12 5/8" x 19 1/8" stainless steel tray and are height adjustable using the Tru-Loc™ friction knob. We offer different base set ups to meet your varying needs.

### Lightweight

**Base:** Two pronged stainless steel base with two hard rubber back casters



### Standard



**Base:** Two pronged stainless steel base with two roller blade back casters



### Four Caster

**Base:** Two pronged stainless steel base with four roller blade casters



Name	Lightweight	 Standard	 Four Caster
Height Adjustment	37 1/4" - 64"	39 1/2" - 62"	39 1/2" - 62"
Model #	1520SS	8841SS	8841MR
Item #	0671520000	0638841000	0638841200
Shipping Weight	26 lbs.	36 lbs.	36 lbs.




# Mayo Stands Stainless Steel

## BENJAMIN MAYO STAND

**Description:** Our Benjamin Mayo stand design is made for ease and durability. Our state of the art double post design minimizes lateral movement and increases overall stability. We offer this design in three different tray sizes.



Via a foot operated mechanism, you can adjust the height from 39 1/2" to 62" high off the ground. A guard over the foot prevents accidental release.

Name	Benjamin Mayo Stands	
<b>Tray Size</b>	<b>12 5/8" x 19 1/8" (Standard)</b>	
Model #	8867SS	
Item #	0668867000	
Shipping Weight	40 lbs.	
<b>Tray Size</b>	<b>16" x 21"</b>	
Model #	8869SS	
Item #	0668869000	
Shipping Weight	41 lbs.	
<b>Tray Size</b>	<b>20" x 25"</b>	
Model #	8871SS	
Item #	0668871000	
Shipping Weight	43 lbs.	

# Mayo Stands

## Construction Specifications

### Construction Specifications

Universal Dimensions	<b>Standard Tray Size:</b> Width: 12 5/8" x Depth: 19 1/8"
Materials	304 stainless steel Chrome
Gauge	All contain stainless steel trays: 20GA <b>Benjamin Mayo Stands:</b> 3/4" diameter, 18GA tubes, 1/8" base <b>Newark Mayo Stands:</b> 1520SS: 1 1/4" diameter, 18GA tube, weighted base 8841SS/8848SS: 1" diameter, 16GA tube, 3/8" base
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	2" plastic hooded twin wheel, hard rubber, or roller blade casters
Notes	<ul style="list-style-type: none"><li>MR conditional Mayo stands available where indicated</li></ul>

### Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at [info@blickman.com](mailto:info@blickman.com).

# **BLICKMAN**

## **OUR PROMISE**

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality. We pride ourselves on providing high-end and personalized service to ensure 100% customer satisfied. If you are not completely satisfied for any reason, we will make it right. All of our equipment is backed by our Limited Lifetime Warranty, so you can be certain that you have made the very best choice for your healthcare facility.

Call Blickman:  
Toll Free: (800) 247-5070  
Local: (973) 330-0557

[www.blickman.com](http://www.blickman.com)  
[info@blickman.com](mailto:info@blickman.com)  
500 U.S. Highway 46 East, Clifton, NJ 07011