

# BLICKMAN



CERVICAL CERCLAGE SET  
CYCLE # 653  
MACHINE NAME: sterilis II  
EXP. date: 11/7 / 2019

L.S.C.S \_ SET  
CYCLE # 653  
MACHINE NAME: sterilis II  
EXP. date: 11/7 / 2019

L.S.C.S-SET#3  
CYCLE # 653  
MACHINE NAME: sterilis II  
EXP. date: 11/7 / 2019

## Case Carts

*Compact, stable, secure transport*

### **Includes:**

- Closed Case Carts
- Open Case Carts
- Shelving Options
- Accessories



# OUR MISSION

## **Customer Focus**

Exceed customer expectations through a commitment to continuous improvement of our products and services.

## **Quality Products**

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

## **Organizational Excellence**

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

## **Marketplace Recognition**

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

## **Collaborative Workplace**


Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.

**BLICKMAN**

# THE BLICKMAN DIFFERENCE

- Operating room surgical staff can transport fully loaded case carts throughout the facility with the doors closed for increased protection and cart cleaning is a simple and easy task.
- We offer options that meet a variety of different surgery needs including neurosurgery, general surgery, outpatient surgery, labor and delivery, and more. We also offer space saving solutions that take up less floor space while being able to handle more extensive procedures.
- Our broad offering of case carts remains the most durable and highest functioning on the market. We are the only manufacturer in our class to offer a lifetime limited warranty. Our advanced engineering capabilities enable us to customize case cart designs to satisfy virtually any requirement your medical facility may have.
- The doors of the closed case carts swing a full 270° and a positive spring acting clip at the cart ends capture doors to stay in the open position during cart washing and while in use in the operating room.

## Key Product Features

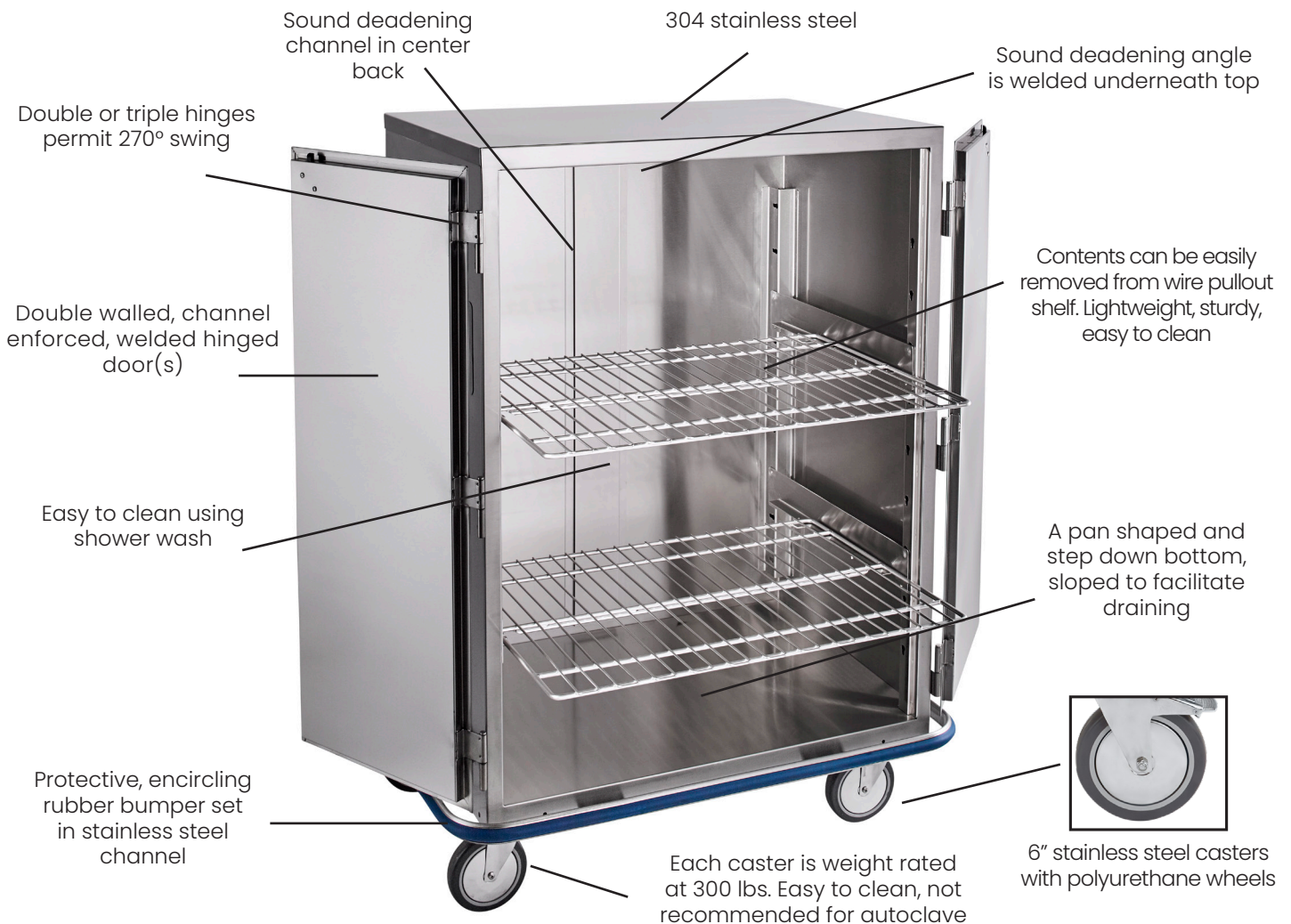

# Case Carts

Maxi to Mini

**Description:** Our white maxi, multipurpose, and mini case carts differ in width, they are all 40 1/2" tall and 29" deep.. They each come standard with one wire shelf. The maxi case cart comes standard with a heavy duty roll-out shelf, and the multipurpose and mini case carts come with a lightweight pull-out shelf.

## Maxi Cart

52"W

Ideal for procedures that require a larger volume of supplies such as neurosurgery and orthopedic surgery

**Internal Dimensions:**  
42 7/8"W x 29 3/8"H x 25"D

## Multipurpose Cart

42" W

Ideal size for most routine surgical procedures, whether or not your facility uses a container system

**Internal Dimensions:**  
32 7/8"W x 29 5/8"H x 25"D

## Mini Cart

29 5/8" W

Designed for labor and delivery and outpatient surgery, where a smaller cart is preferred

**Internal Dimensions:**  
20 1/2"W x 29 5/8"H x 25"D

### Push Handle



Name

Maxi

Multipurpose

Mini

### Solid Doors

	Maxi	Multipurpose	Mini
Model #	CCC1-19	CCC2-19	CCC3-19
Item #	2293331000	2293332000	2293333000
Floor Space	10 1/2 ft <sup>2</sup>	8 1/2 ft <sup>2</sup>	6 ft <sup>2</sup>
Shipping Weight	266 lbs.	225 lbs.	168 lbs.

### Glass Doors

	Maxi	Multipurpose	Mini
Model #	CCC1-19G	CCC2-19G	CCC3-19G
Item #	22933G1000	22933G2000	22933G3000
Floor Space	10.47 ft <sup>2</sup>	8.46 ft <sup>2</sup>	5.97 ft <sup>2</sup>
Shipping Weight	260 lbs.	219 lbs.	165 lbs.

**Custom Designs Available Upon Request**

Please email us at [info@blickman.com](mailto:info@blickman.com)

# Case Carts

Space Saving

**Description:** Our space saver case carts are designed to give a greater capacity while using the same floor space as the mini and multipurpose case carts. Their height and additional shelving provide additional storage capacity. Instead of 40 1/2" tall, they are 55 1/2" tall, and instead of one shelf, they come standard with two lightweight pull-out wire shelves.

## Space Saving Cart

29 5/8" W

This cart can handle most routine surgical procedures.

### Internal Dimensions:

20 1/2"W x 44 5/8"H x 25"D



## Ultra Space Saving Cart

42" W

This cart can handle more extensive procedures.

### Internal Dimensions:

32 7/8"W x 44 5/8"H x 25"D



## Recessed Handle



Name	Space Saving	Ultra Space Saving
<b>Solid Doors</b>		
Model #	CCC4-19	CCC5-19
Item #	2293334000	2293335000
Floor Space	6 ft <sup>2</sup>	8 1/2 ft <sup>2</sup>
Shipping Weight	208 lbs.	276 lbs.
<b>Glass Doors</b>		
Model #	CCC4-19G	CCC5-19G
Item #	22933G4000	22933G5000
Floor Space	6 ft <sup>2</sup>	8 1/2 ft <sup>2</sup>
Shipping Weight	198 lbs.	262 lbs.

**Custom Designs Available Upon Request**

Please email us at [info@blickman.com](mailto:info@blickman.com)

# Case Carts

Open

**Description:** Our open case cart designs are easy to clean, dry, load, and maneuver. The full encircling bumper and solid bottom prevents cross contamination. This design is ideal for carts that go directly to a sterile core, bypassing non-sterile transport.

## Overall Dimensions:

Width: 42" x Height: 42" x Depth: 29"



Name	Solid Shelves	Wire Shelves
Model #	OCC3	OCC4
Item #	2410396000	2410496000
Shipping Weight	219 lbs.	141 lbs.

**Custom Designs Available Upon Request**  
Please email us at [info@blickman.com](mailto:info@blickman.com)



# Case Carts

Interior Shelving

## STANDARD

### Wire Roll Out Shelf



CCC1

### Wire Pull Out Shelf



CCC2, 2E, 3, 4, and 5

## OPTIONAL



### Roll Out Shelving

### Pull Out Shelving

Case Carts	Wire	Perforated	Solid	Wire	Perforated	Solid
CCC1	611C1RWKIT	611C1RPKIT	613C1RSKIT	N/A	61517PIKIT	6152042KIT
CCC2/CCC5	610C2RWKIT	610C2RPKIT	610C2RSKIT	6152272KIT	6152275KIT	6152274KIT
CCC3/CCC4	610C1RWKIT	610C1RPKIT	610C1RSKIT	6152273KIT	6152271KIT	6152073KIT

# Case Carts

Exterior Shelving

## EXTENSION SHELVES



Compatible with CCC2 Only

- Expand your work surface area -

## DROP SHELVES



Compatible with CCC1, 2, & 3

- Expand your work surface area -

**Custom Designs Available Upon Request**

Please email us at [info@blickman.com](mailto:info@blickman.com)



# Case Carts

## Accessories

### ACCESSORIES

#### Clean/Dirty Indicators



US Design Patent #: US D883,378 S  
US Utility Patent#: US 10,643,505 B1



US Design Patent #: US D883,379 S  
US Utility Patent#: US 10,643,505 B1

Part #	Description	Shipping Weight
5278628KIT	Angle Model	2 lbs.
5200045KIT	Dial Model	2 lbs.
5278618KIT	Pod Model	2 lbs.



#### Card Holder

We offer a 8.5" x 11" card holder. Custom sizes and placement options available. Stainless steel frame withstands prolonged use and cutouts allow for drainage.



#### Tote Boxes

Sturdy, long lasting tote boxes are ideal for organizing and storing small supplies that may slip through the spaces of wire shelves.

Part #	Description	Shipping Weight
5206701185	Card Holder	N/A - Factory Installed
9854194720	Tote Boxes	5 lbs.



# Case Carts

## Construction Specifications

### Construction Specifications

Universal Dimensions	<b>Closed Case Carts:</b> 29" deep, varying widths. Standard Height: 40 1/2"; Space saving height: 55 1/2" <b>Open Case Carts:</b> 42" x 42" x 29"
Materials	304 stainless steel, non-magnetic
Gauge	<b>Closed Case Carts:</b> Body & Doors: 18GA; Top & Chassis: 16GA <b>Open Case Carts:</b> Top & Bottom: 16GA; Shelf: 18GA <b>Bumper Channel:</b> 14GA
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	Four 6" stainless steel swivel casters with polyurethane wheels, weight rated at 300lbs (one with directional lock)
Notes	Found on next page

### Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at [info@blickman.com](mailto:info@blickman.com).

**Custom Designs Available Upon Request**  
Please email us at [info@blickman.com](mailto:info@blickman.com)

# Case Carts

## Construction Specifications

---

### Notes Continued

---

#### **CLOSED CASE CARTS:**

- **Body:** One piece top and one piece body. Sound deadening angle welded under the top. Sound deadening channel on the center back. A pan shaped, step down bottom with a slope that facilitates draining. Horizontal intersections filled with heat resistant silicone. The body shell is bolted to the chassis. A tubular push handle is provided at one end of the cart.

- **Doors:** Double walled, channel enforced, (double or triple) hinged doors with paddle latch pulls on leading edges. Doors less than 36" high swing on two hinges, while doors greater than 36" high swing on three hinges. All doors permit full 270° swing. In fully opened position, positive spring acting clips at the cart ends capture doors for cart washing or operating room use.

- **Shelves:** Shelves slide on adjustable stainless steel wires, suspended from keyhole punched standards welded to the cart interior.

#### **OPEN CASE CARTS:**

- **Frame:** U-shaped tubular framing formed around 3" radii or tubular corner posts stabilized and capped by a pan shaped top welded to posts at contact points. Posts are fitted with casters and either rolling bumpers or encircling rubber bumpers.

- **Tops, bottoms, fixed shelves:** Fabricated with turned down edges, either hemmed or flanged, with corners notched to sit over or snugly against the inside contour of the upright members. Notched corners are fully welded to the uprights the full height of the facing edges as well as around the notch.

- **Adjustable shelves:** Solid rods welded between the tubular uprights at each end, stabilizing the framing members while providing varying levels of shelf adjustment.

---

**Custom Designs Available Upon Request**

Please email us at [info@blickman.com](mailto:info@blickman.com)

# **BLICKMAN**

## **OUR PROMISE**

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality. We pride ourselves on providing high-end and personalized service to ensure 100% customer satisfied. If you are not completely satisfied for any reason, we will make it right. All of our equipment is backed by our Limited Lifetime Warranty, so you can be certain that you have made the very best choice for your healthcare facility.

Call Blickman:  
Toll Free: (800) 247-5070  
Local: (973) 330-0557

[www.blickman.com](http://www.blickman.com)  
[info@blickman.com](mailto:info@blickman.com)  
500 U.S. Highway 46 East, Clifton, NJ 07011