

SmartGown™ Surgical Gowns

All around protection you need. Comfort you want.



SmartGown is the No. 1 brand¹ of AAMI Level 4 surgical gowns in the U.S.

Our SmartGown Surgical Gowns feature complete AAMI Level 4 protection with a unique RESPONSIVEMEMBRANE™ technology to transport moisture from the skin to the surface of the fabric — keeping clinicians cool, dry and comfortable

SmartGown™ Breathable Surgical Gowns			
Sterile	BNS	Description	Qty. Cs
89005	89005N	Small/medium	20
89015	89015NA	Large	20
89045	89045NA	X-large	18
89075	89075NA	XX-large	14
89085	89085N	XXX-large	12
89095		XXXX-large	12
SmartGown™ Breathable Surgical Gowns with raglan sleeves			
39015	39015NA	Large	20
39019	39019N	Large, X-long	14
39045	39045NA	X-large	18
39049	39049NA	X-large, X-long	14
39075	39075NB	XX-large	14
39079	39079NB	XX-large, X-long	14
39099	39099NA	XXX-large, X-long	14
32474	32474NA	X-large, X-long, A-line	16
SmartGown™ Breathable Scrub Nurse* Gowns			
	324740NA	X-large, X-long, A-line	n/a
	890015NA	Large	n/a
	890045NA	X-large	n/a

AAMI Levels determined per ANSI/AAMI PB70:2012

¹ 2016 US GHX sales and units data across all GHX markets

* Scrub Nurse Gowns are bulk non-sterile (BNS) gowns intended for procedure kits. The gown and hand towel are wrapped in SMS material for placement in the kit.



Material attributes:

SFS gown base fabric

SFS: SMS (Spunbond/Meltblown/Spunbond) /Polyester film/SMS.

Spunbond and meltblown composition is polypropylene.


Sleeve

SFS (Spunbond/Polyester film/Spunbond).

Spunbond composition is polypropylene.

White cuff

Polyester

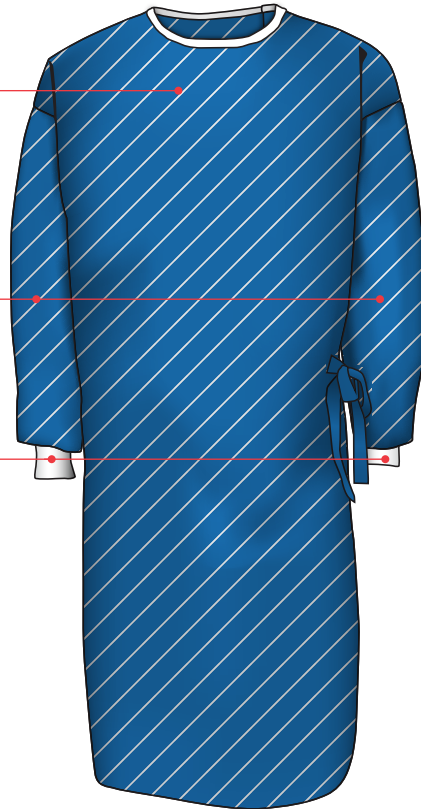

Drawing is for illustrative purposes only.

Biocompatibility:

Cardinal Health™ Surgical Gowns are tested for and meet the toxicity requirements contained in the ISO 10993 standards for biological evaluation and biocompatibility testing of medical devices.

Flammability standard:

Cardinal Health™ Surgical Gown materials meet the "Standard for the Flammability of Clothing Textiles" for disposable surgical drapes and gowns required by the Consumer Product Safety Commission (16 CFR Part 1610, formerly CS 191-53). The base fabric comprising Cardinal Health™ Surgical Gowns meets Class 1: Normal Flammability, which is the rating required for consumer use. Cardinal Health promotes the safe use of our products and does not recommend placing any ignition or heat source in contact with our surgical gowns.

