

model **6500**  
**WHEELCHAIR  
SCALES**

**TWO-WAY EASY-ACCESS RAMPS**



- Robust 1,000 lb / 450 kg capacity
- Low-profile platform for easy patient access
- Integral dual wheels and carrying handles for transport
- EMR/EHR-ready serial, USB, and optional Ethernet connectivity

Detecto  
Model 6500  
Low-Profile Portable  
Wheelchair Scale



**With built-in ramps for entering and exiting the scale,** Detecto's 6500 portable wheelchair is perfect for clinical weighing applications. For added convenience, you may enter the wheelchair weight via pushbutton or keypad tare, so that only the patient weight is displayed. With a high capacity of 1,000 lb / 450 kg, this scale is ideal for bariatric weighing needs.

**Versatile Weighing Options with the 6500:**

Wheelchairs • Standing • Straight-legged Chairs • Walkers



- 1 Easy Transport**
- 2 Lower into Position**
- 3 Ready for Weighing**

The 6500's indicator secures firmly to the base for safer movement and storage.



**Carrying Handles** ease transport and allow the scale to be used in multiple locations.



**Integral Wheels** combined with a lightweight frame and battery power make the 6500 completely mobile.

**Optimum patient weighing "on-the-move"**

## EFFICIENT MEASUREMENT



The 6500's digital weight indicator is the clinical industry's most advanced, featuring EMR/EHR-ready serial, USB and optional Ethernet interface, BMI calculation, yellow "quick keys" for basic weighing, and a multi-line LCD display. The scale is powered by 6 C batteries (not included) or optional AC adapter.



Detecto's 6500 scale features patient ID barcode scanning capability which sends a 9-digit patient ID to the indicator for efficient EMR/EHR recordkeeping.


The weight indicator features 19 keys, 5,000 divisions, battery or AC power, and selectable key lockout.

- Up to 99 tares (optional) to store wheelchair weights
- Hold/Release key to retain weight after leaving the scale
- Non-skid rubber mat for patient assurance and safety




**USA-Made Clinical Quality Construction**

model **6500**  
**WHEELCHAIR**  
**SCALES**

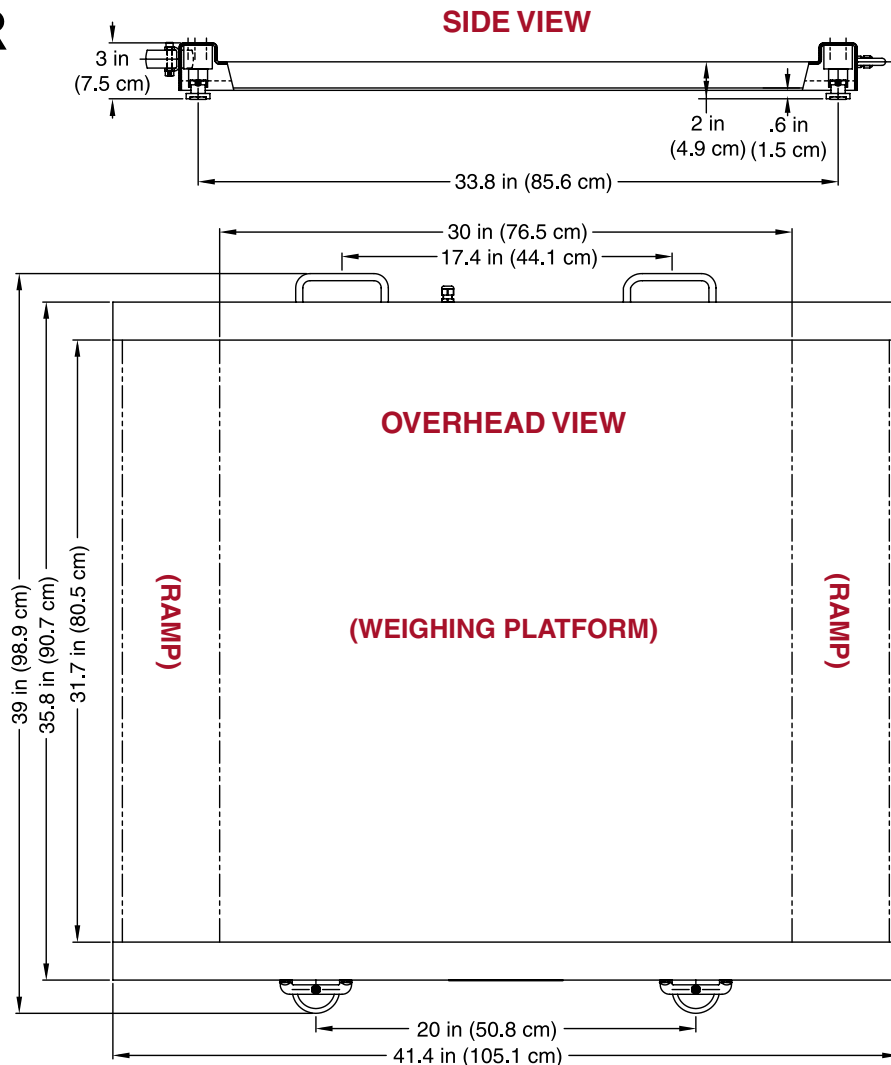

<b>Capacity</b>	1,000 lb x 0.2 lb / 450 kg x 0.1 kg
<b>Display</b>	Weight: Five digit, seven segment .5 in / 13 mm H LCD Height: Four digit, seven segment .375 in / 10 mm H LCD
<b>Keyboard</b>	Membrane type with 19 keys
<b>Connectivity Ports</b>	RS232 serial, USB, and optional Ethernet
<b>Display Dimensions</b>	9.5 W x 6.6 H x 2.1 D in / 24 W x 17 H x 5 D cm
<b>Power</b>	6 C cell Ni-Cad or NiMH batteries (not included) OR an optional 100 to 240 VAC 50/60Hz 12 VDC 1A wall plug-in UL/CSA listed AC power adapter (Detecto part number 6800-1045)
<b>Resolution</b>	5,000 divisions, maximum
<b>Operating Temperature</b>	14° to 104° F (-10° to +40° C)
<b>Enclosure Construction</b>	Aluminum
<b>UPC Code</b>	809161139204
<b>Shipping Weight</b>	65 lb / 29.4 kg

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

 **DETECTO**<sup>®</sup>  
 A Division of Cardinal Scale Manufacturing Co.  
 203 E. Daugherty, Webb City, MO 64870 USA  
 Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001  
 www.Detecto.com

**Sold By:**

