

Model DS1125

Infant Measuring Rod



The DS1125 infant measuring rod provides an accurate measurement tool that can be moved to the patient or mounted on a wall or pediatric examination table. Precision engineered and constructed of heavy gauge aluminum, the Doran Model DS1125 infant measuring rod will deliver accurate height readings for years and years. Simply flip the measurement arms away from the measuring rod and slide the measurement arm into position. The patient's length is highlighted with an easy to read, bright red line.

Marked in millimeters and inches, the Model DS1125 infant measuring rod can be used to measure adult or pediatric patients.

MEASURING RANGE

0 - 39" (0 - 99 cm)

GRADUATIONS

1/16" (1 mm)

CONSTRUCTION

Aluminum and ABS Plastic

HEIGHT ROD DIMENSIONS


40.5" L x 1.5" W x 0.5" H

MEASUREMENT ARM DIMENSIONS

10" L x 1.5" W x 1.5" H

WARRANTY

2 Two Years

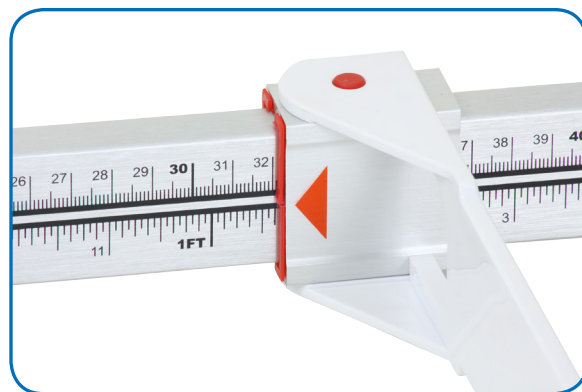

The precise markings of the DS1125 ensure accurate infant length data.



Doran Scales, Inc. 1-800-264-4106
1315 Paramount Parkway Fax: 630-879-0073
Batavia, IL 60510 www.doranmedical.com