

# 25

S E R I E S

A New Standard for Portable Therapeutic Ultrasound



DYNATRON®  
125



Touch Screen



Plug N' Play



Battery Option



Portable



Slim Input



3-Frequency



Modifiable




3 Head Option






## Touch Screen User Interface

The modern innovation of the touch screen makes the Dynatron 125 the ultimate device in Ultrasound treatment delivery. The user experience is maximized as a result of the simple and easy to use interface. With the touch of a button, an Ultrasound treatment is deliverable in seconds.

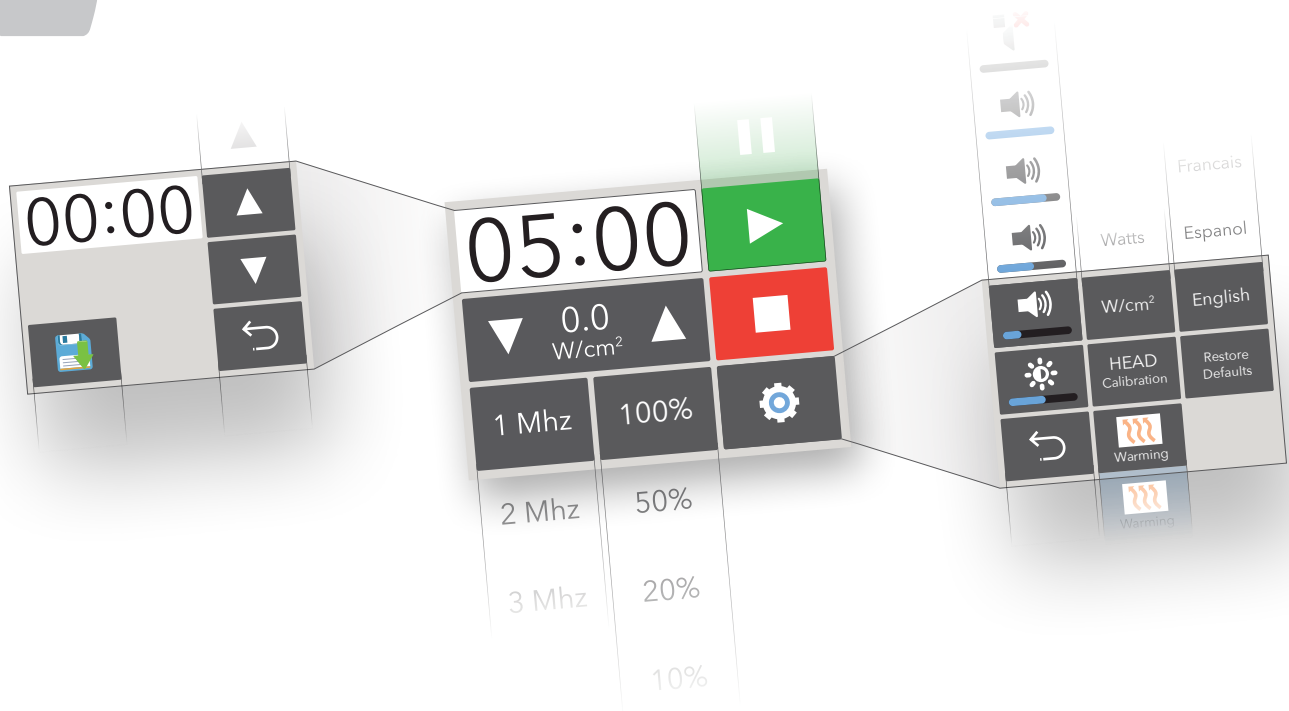

## Electrical Stim Input

This jack is used for combining electrical stimulation with Ultrasound. The output end of a lead connected to a Dynatron electrotherapy device may be plugged into this jack. The current from the stim device is delivered through the soundhead.

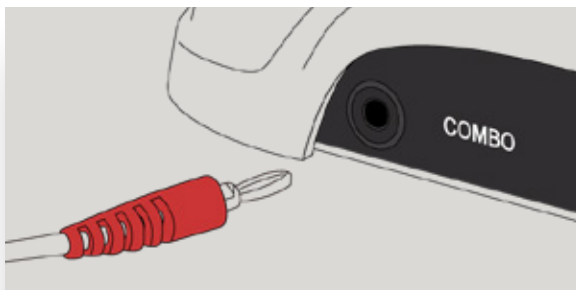

## Optional Battery

The Dynatron 125 may be operated with power supplied by an external 12-volt battery (not included), offering greater flexibility and portability in treatment sites.





## Three Optional Soundheads

Dynatronics offers Ultrasound SmartHeads in three sizes, 2 cm<sup>2</sup>, 5 cm<sup>2</sup>, and 10 cm<sup>2</sup>. One soundhead of your choice is included with the purchase of a Dynatron 125.



## SmartHead Technology

The Dynatron 125 features SmartHead technology, having all calibration information stored within the soundhead. Whichever soundhead is plugged into the device will be uniquely recognized and be ready to use.



## 3-Frequency Ultrasound

With the Dynatron 125 you can select three different frequencies without changing the soundhead. The 2 cm<sup>2</sup>, 5 cm<sup>2</sup>, and 10 cm<sup>2</sup> soundheads each operate at 1, 2, and 3 MHz, providing the greatest flexibility in depth of treatment.



WSH02  
2 cm<sup>2</sup> Ultrasound SmartHead

WSH05  
5 cm<sup>2</sup> Ultrasound SmartHead

WSH10  
10 cm<sup>2</sup> Ultrasound SmartHead



## Portable Modalities

The small light-weight design of the Dynatron 125 makes it easy to carry and move around the clinic. The battery powered option allows the device to travel to any location, providing Ultrasound treatments immediately, wherever they are needed.



## In-Treatment Modifiable

Treatment time, power, frequency, or duty cycle can be modified while a treatment is in progress—without stopping the treatment.



## Additional Features

Customizable Duty Cycle from 10, 20, 50 percent, and Continuous.  
Head warming (optional selection). When the device is ON but idle, the soundhead is warmed automatically.  
Two-Year Warranty.

25  
S E R I E S



DYNATRON®  
125

D125B  
Dynatron 125 Portable Ultrasound



## Technical Specifications

Meets IEC 60601-1 for Safety and IEC 60601-2-5 for Ultrasound Performance.

Supply Voltage: 110-240 Vac  
Power: 65 watts  
Environment Temperature: +10° C - 30° C  
Humidity: 30% - 50%

Fuse: 1.6A  
Dimensions: 9" x 7" x 4"  
Weight: 4 lbs.



7030 Park Centre Drive | Salt Lake City, UT | 84121  
toll-free [800] 874-6251 | local [801] 568-7000 | fax [801] 568-7711  
[www.dynatronics.com](http://www.dynatronics.com) | email: [sales@dynatronics.com](mailto:sales@dynatronics.com)