

# 25

S E R I E S

Electrotherapy | Ultrasound | Combination



Plug N' Play



Battery Option



3 Head Option



Stim Input



3-Frequency



Modifiable



Target





# Dynatron® 25 Series




## 3 & 5 Channel Models with Unequaled ComboStim™ Combination Therapy

The only 3 and 5 channel devices available. Channels allow fully-independent treatment setups offering Interferential, Premodulated, High Volt, Biphasic, Russian and Fixed Frequency IFC/Premod. In addition, the 925 and 825 models include Dynatronics' Ultrasound Comboplus feature with the power to deliver up to 5 channels of stim and ultrasound – all at the same time.



## Target & Target Sweep for Interferential Treatments

It is easy and efficient to focus treatment precisely where it's needed! Simply glide your finger across the TARGET touchpad to move the center of interference to the site of your patient's pain. Available on all 25 Series electrotherapy devices.

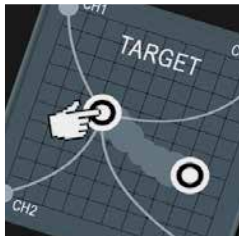

**WSH10**  
10 cm² Ultrasound SmartHead

**WSH05**  
5 cm² Ultrasound SmartHead

**WSH02**  
2 cm² Ultrasound SmartHead



## 3 Frequency Ultrasound

Dynatron 25 Series models offer 1, 2, and 3 MHz frequencies providing great flexibility for depth of treatment. Choose 1 MHz for deep treatments, 2 MHz for moderate depth, or 3 MHz for superficial depth.



## 3 Optional Soundheads

Dynatronics offers Ultrasound SmartHeads in three sizes, 2 cm², 5 cm², and 10 cm². One 5 cm² soundhead is included with the purchase of any 925, 825, or 125 model.



## SmartHead Technology

The 25 Series features Ultrasound SmartHead technology, that stores calibration information within the soundhead. Whichever soundhead is plugged into the device will be uniquely recognized and ready to use.



## In-Treatment Modifiable

Treatment time, power, frequency, or duty cycle can be modified while a treatment is in progress—without stopping the treatment.



	<b>Dynatron 925</b> (D925T)	<b>Dynatron 825</b> (D825T)	<b>Dynatron 625</b> (D625T)	<b>Dynatron 525</b> (D525T)	<b>Dynatron 125</b> (D125B)
IFC	✓	✓	✓	✓	
Premod	✓	✓	✓	✓	
Fixed Frequency*	✓	✓	✓	✓	
Biphasic	✓	✓	✓	✓	
Russian	✓	✓	✓	✓	
High Volt	✓	✓	✓	✓	
Combo**	✓	✓			
Ultrasound	✓	✓			✓
Channels	5	3	5	3	Stim Input Jack
Battery Option	✓	✓	✓	✓	✓
<b>Technical Specifications</b>					
Supply Voltage	100-240 V-, 50/60 Hz				110-240 Vac
Power	100 watts				65 watts
Environment Temp.	10° C - 30° C				10° C - 30° C
Humidity	30% - 75%				30% - 50%
Fuse	250 V, T5AH High Break				1.6A
Dimensions	18.5" W (41.91cm) x 4.0" H (10.2cm) x 12.0" D (30.5cm)				9" x 7" x 4"
Weight	8.2 lbs.				4 lbs.

\* IFC/Premod \*\* Electrotherapy / Ultrasound

