


## THE 3001KL-AM

### ANTIMICROBIAL DIGITAL PLATFORM SCALE

The all-digital 3001KL-AM supports up to 1000 pounds and provides the optimum quality, reliability, and ease of use with a state-of-the-art digital height rod. Designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities, The scale features a large non-skid low profile platform, "live" handrails for patient stability, and rugged scuff-proof wheels that enable portability. The scale also offers an antimicrobial keypad, handrails, and transport handle to give patients and providers additional safety. The scale's high-contrast color screen provides a user-friendly interface that includes BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently

#### Antimicrobial

Keypad  
Handrails  
Transport Handle

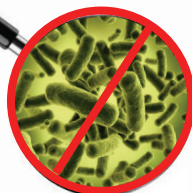

The 3001KL-AM offers two optional height rod options with quick and easy to assemble brackets, a digital height rod 245EHR-3001 or a mechanical height rod 3001HR. The 245EHR-3001 market-leading digital height rod provides the most accurate height measurement in its class, with a graduation/resolution of 1mm/0.1" and a full measuring range of 11 3/4" – 79 1/2" / 30cm – 202cm and a digital range of 3ft 11.2" – 6ft 7.5" / 120 cm - 202 cm. The 3001HR features a durable height rod construction and an easy-slide head-piece that allows caregivers to quickly and effortlessly measure and read a patient's height. The 3001HR provides a measuring range of 2"-90" in 1/16" increments (6 cm-230cm in 1 mm increments).

The scale provides easy integration and seamless communication with today's EMR applications. The 3001KL-AM includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

### THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 3001KL-AM is your choice for safety and reliability.

**120+** LOCAL  
FIELD SPECIALISTS

**100%**  
CONNECTIVITY

**#1** TRUSTED  
BRAND

**100 YEARS**  
OF SERVICE



**Health o meter®**  
Professional

Health o meter® Professional 3001KL-AM

weigh**easier**®

## TECHNICAL SPECIFICATIONS

### Capacity

1000 lb / 454 kg

### Resolution

0.2 lb / 0.1 kg

### EMR Connectivity

Via USB, optional Pelstar® wireless technology

### Display

1 ½" / 38 mm High-Contrast Color TFT-LCD Screen

### Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

### Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare  
Pre-tare  
Hold/release  
Reweigh  
Recall  
Auto zero  
Auto off  
Time/date  
Variable auto off time  
Audible/mute sound option  
Parent/child weighing function  
Everlock® - permanent weight lock

### Platform Size (w x d x h)

22" x 15 ¾" x 2 ⅝"  
(559 mm x 400 mm x 64 mm)

### Product Footprint (w x d x h)

22" x 33 ⅛" x 50"  
(559 mm x 841 mm x 1270 mm)

### Product Weight

46 lb / 21 kg

### Portability

2 wheels for ease of movement

### Warranty

2-year limited manufacturer warranty

## ORDERING INFORMATION

MODEL	DESCRIPTION
3001KL-AM	Antimicrobial digital platform scale, kilograms and pounds
3001KG-AM	Antimicrobial digital platform scale, kilograms only
3001KL-AM-BT	3001KL-AM with built-in Pelstar wireless technology
3001KG-AM-BT	3001KG-AM with built-in Pelstar wireless technology
3001KL-AMHR	3001KL-AM with mechanical height rod
3001KG-AMHR	3001KG-AM with mechanical height rod
3001KL-AMHR-BT	3001KL-AMHR with built-in Pelstar wireless technology
245EHR-3001	Optional digital height rod for 3001 models*
3001HR	Optional mechanical height rod for 3001 models
SS-3001KL-AM	ScaleSurance 2 year extended warranty (3001 models)**
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
PELSTARONE-P	Wireless thermal printer

\* Compatible with 3001 series scales including a date code with the last two digits of 21 or higher. An upgrade option is available for models including a date code with the last two digits of 17 or higher. Not available for AMX series scales.

For assistance, contact Customer Service at 800-815-6615.

\*\* Not available in all countries

### Optional Height Rods Available



To learn more, call **1-800-815-6615**  
or visit **homscales.com/3001KL-AM**

**Health o meter®**  
Professional




# Antimicrobial Digital Platform Scale with Extended Handrails

3001KL-AMX

antimicrobial  
**Health o meter®**  
Professional  
scales

antimicrobial technology

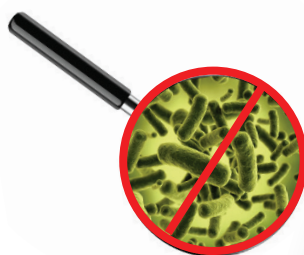

Patented\* Antimicrobial Powder Coating uses sustainable silver ions to inhibit the growth and transmission of harmful microbes on treated surfaces.

\*Patents US6093407(A), US643241 (B1)

## FALL PREVENTION

***Safe patient handling and reducing costs associated with inpatient falls are primary concerns for healthcare providers***

### ✔ **“Live” Handrails**

Advanced technology allows the patient to hold onto handrails for stability while the scale accurately calculates weight

### ✔ **Extended Handrails**

Longer length handrails provide extra support and help patients maintain balance when stepping onto the platform

### ✔ **Low Profile Platform**

requires only minimal step up, offering the patient a more comfortable and stable weighing experience

## WORKFLOW EFFICIENCY

***Simple to use and intuitive scales allow caregivers to focus on patients instead of complicated equipment***

### ✔ **On-Screen Help Menu**

Guides users through operating functions and provides easy access to multiple settings

### ✔ **Interactive Display**

High contrast color screen provides a user-friendly interface

### ✔ **EMR Capable**

Scale has the ability to interface with a PC or other electronic device via USB or optional Pelstar wireless technology

### ✔ **Everlock® (patent pending)**

Option that allows users to permanently lock the scale in either LB or KG to improve safety



**LIVE HANDRAILS  
FOR PATIENT STABILITY**

## INFECTION PREVENTION

***Reducing patient's risk against Health-care Associated Infections is a major challenge for medical facilities***

### ✔ **Antimicrobial\* Powder Coating**

The silver ionic technology suppresses the growth of bacteria and helps protect against cross-contamination

### ✔ **Antimicrobial Keypad**

The keypad features antimicrobial protection

## DURABILITY

***Facilities require a product to withstand frequent use, movement, and disinfection***

### ✔ **Convenient Transport Handle**


allows the user to easily move the scale to where it is needed in the facility

### ✔ **Rugged Scuff-Proof Wheels**

enable smooth and easy portability

# Antimicrobial Digital Platform Scale with Extended Handrails

3001KL-AMX

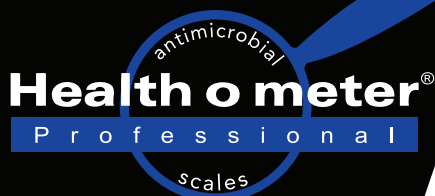

antimicrobial technology

CONVENIENT  
TRANSPORT  
HANDLE



## SPECIFICATIONS

- Capacity: 1,000 lb / 454 kg
- Resolution: 0.2 lb / 0.1 kg
- Platform Size (w x h x d):  
22" x 15 3/4" x 2 5/8" / 559 mm x 400 mm x 67 mm
- Handrail Width: 20" / 508 mm
- Product Footprint (w x h x d):  
22" x 34 1/2" X 52" / 559 mm x 876 mm x 1321 mm
- Product Weight: 43 lb / 19.5 kg
- Display: 1 1/2" / 38 mm High Contrast Color TFT-LCD Screen
- Rugged Scuff-Proof Wheels: 2
- Power Source: 100-240V Adapter ADPT30 (included), 6 D batteries (not included)
- Connectivity: USB, optional Pelstar wireless technology
- Functions: LB/KG Conversion, LB/KG Lock Out, Body Mass Index (BMI), Zero, Tare, Pre-tare, Hold/Release, Recall, Reweigh, Auto Zero, Auto Off, Time/Date, Variable Auto Off Time, Audible/Mute Sound Option, Everlock®
- 2 year limited warranty

Model	Description
3001KL-AMX	Digital platform scale with extended handrails, pounds & kilograms
3001KG-AMX	Digital platform scale with extended handrails, kilograms only
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
SS-3001KL-AMX	ScaleSurance 2 year extended warranty for 3001KL-AMX & 3001KG-AMX models*
PELSTARONE-P	Wireless thermal printer

