

LOOK™ SUTURE

PRODUCT CATALOG



- Consistency
- Quality
- Value

WELCOME

The LOOK™ line of sutures is one of the best known and most trusted suture brands on the market. From the beginning, these sutures have been known for their quality—and that legacy of quality continues to meet the demands of surgery professionals worldwide. The precision manufacturing of our LOOK™ brand sutures ensures consistent needle strength and durability that holds up to every pass.

LOOK™ Sutures are offered in a variety of absorbable and non-absorbable materials and are attached to needles with superior sharpness, penetration, and strength characteristics.

Surgical Specialties Corporation is proud to provide the LOOK™ suture, along with flexible and accommodating customer service, and value to its customers.

Table Of Contents

How to Use Our Catalog	1
Materials Chart	2
A Guide To Surgical Sutures	3
Needle Anatomy and Terminology	4
Needle Shapes And General Uses	5
Needle Points and Cutting Edges	6
Needle Charts	7
Absorbable Sutures	
Plain Gut.....	9
Chromic Gut	11
PolySyn FA™ (polyglycolic acid)	13
PolySyn™ (polyglycolic acid)	14
Non-Absorbable Sutures	
Silk	18
Nylon.....	21
Polypropylene.....	24
Polyviolene™ (Polyester) - Uncoated.....	26
Additional Products	27
Product Index.....	31

Introduction

How to Use Our Catalog

The pages in this catalog show the following information for each product code. Some products feature an expanded chart with detailed measurements (not shown in this sample) to ensure proper needle selection.

Needle Description - Contains an illustration of the needle. The LOOK™ needle code is purple, the comparable Ethicon™ needle code is black, and the comparable Syneture™ needle code is dark green.


Needle Geometry - Type of needle with tip geometry

Material - Color and style of suture

Precision
Reverse Cutting



Length - Length in inches/centimeters

Needle	Material	Length	Size				
			5-0	4-0	3-0	2-0	0
PC31 19mm 3/8 Circle  PS-2 P-12	Undyed Braid	18"/45cm	495B J495G SL5636G	496B • SL5627G			

Size Grid - Indicates the U.S.P. size for each product code listed.

• • indicates that other manufacturers do not have a comparable product code or needle. LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

LOOK™ product code

Ethicon™ comparable product code




Syneture™ comparable product code

Materials Chart






LOOK™ Suture	Ethicon™	Syneture™	Material
Plain Gut	Surgical Gut	Plain Gut	Natural Absorbable Gut
Chromic Gut	Surgical Gut	Chromic Gut	Natural Absorbable Gut
PolySyn FA™	Vicryl Rapide™	Velosorb™ Fast	Synthetic Absorbable, Braided, Coated Polyglycolic Acid
PolySyn™	Vicryl™	Polysorb™/Dexon™	Synthetic Absorbable, Braided, Coated Polyglycolic Acid
Silk	Permahand Silk™	SofSilk™	Non-Absorbable, Braided, Coated Silk
Nylon	Ethilon™	Monosof™/Dermalon™	Non-Absorbable, Monofilament Nylon
Polypropylene	Prolene™	Surgipro™ II	Non-Absorbable, Monofilament Polypropylene
Uncoated Polyviolene™	Mersilene™	Surgidac™	Non-Absorbable, Braided, Uncoated Polyester
Coated Polyviolene™	Ethibond™	Ti-Cron™	Non-Absorbable, Braided, Coated Polyester

A Guide To Surgical Sutures

Absorbable


Material	Material Description	BSR Profile	Mass Absorption
 <p>• Surgical Gut Suture - Plain/Chromic * Plain Gut/Chromic Gut</p>	<p>Twisted</p> <ul style="list-style-type: none"> Plain & Chromic Packaged Dry or in Wetting Solution <p>Common uses: General Closure, OB/Gyn, Urology, Bowel Anastomosis</p>	<p>In vivo strength retention: Varies</p>	<p>Varies depending upon type of suture (plain vs. chromic) and site of placement</p>
 <p>• VICRYL RAPIDE™ * Velosorb™ Fast</p>	<p>Braided</p> <ul style="list-style-type: none"> Undyed (clear) • Coated <p>Common uses: General Closure, OB/Gyn, Urology</p>	<p>In vivo strength retention: 48% at 7 days 0% at 14 days</p>	<p>Essentially complete between 49 and 63 days</p>
 <p>• VICRYL™ * Polysorb™ / Dexon™</p>	<p>Braided</p> <ul style="list-style-type: none"> Undyed (clear) • Dyed (violet) <p>Monofilament</p> <ul style="list-style-type: none"> Undyed (clear) • Dyed (violet) Coated & Uncoated <p>Common uses: Skin Closure, OB/Gyn, Urology, Gastrointestinal</p>	<p>In vivo strength retention: 50% at 14 days 20% at 21 days</p>	<p>Essentially complete between 50 and 90 days</p>

Non-Absorbable

Material	Material Description	Qualities
 <p>• PERMA-HAND™ Silk * SofSilk™</p>	<p>Braided</p> <ul style="list-style-type: none"> Undyed (white) • Dyed (black) • Coated <p>Common uses: General & Skin Closure, Gastrointestinal, Cardiovascular Surgery, Plastic Surgery, Ophthalmology</p>	<ul style="list-style-type: none"> Excellent handling and tying characteristics
 <p>• ETHILON™ * Monosof™/Dermalon™</p>	<p>Monofilament</p> <ul style="list-style-type: none"> Undyed (clear) • Dyed (blue) (black) Packaged dry or in wetting solution <p>Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery, Gastrointestinal</p>	<ul style="list-style-type: none"> Passes easily through tissue
 <p>• PROLENE™ * Surgipro™/Surgipro™ II</p>	<p>Monofilament</p> <ul style="list-style-type: none"> Undyed (clear) • Dyed (blue) <p>Common uses: General & Skin Closure, Neurosurgery, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery</p>	<ul style="list-style-type: none"> Inert Easy to handle Passes easily through tissue
 <p>• MERSILENE™ * Surgidac™</p>	<p>Braided</p> <ul style="list-style-type: none"> Undyed (white) • Dyed (green) Uncoated <p>Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures</p>	<ul style="list-style-type: none"> Smooth Easy to handle Good knot security
 <p>• ETHIBOND EXCEL™ * Ti-Cron™</p>	<p>Braided</p> <ul style="list-style-type: none"> Undyed (white) • Dyed (green) Coated <p>Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures</p>	<ul style="list-style-type: none"> Smooth Easy to handle Good knot security

• Comparable Ethicon™ Material * Comparable Synture™ Material

Needle Anatomy And Terminology

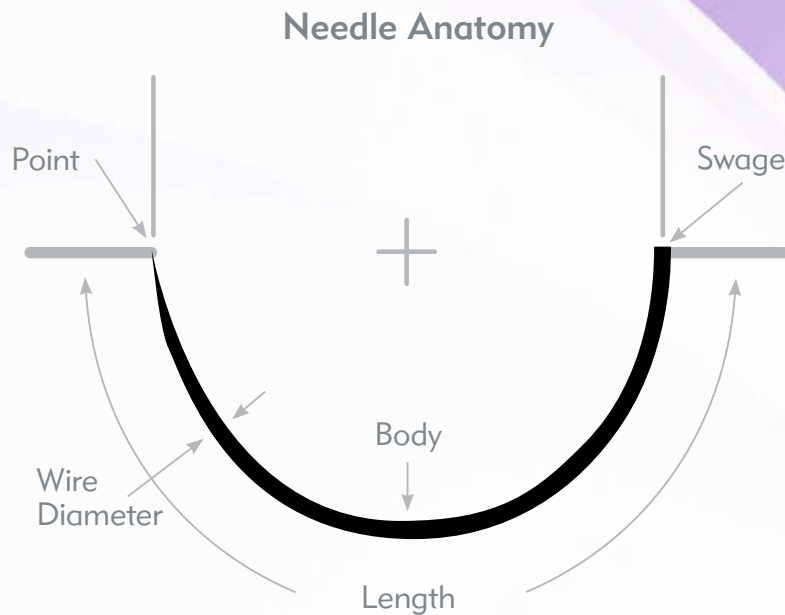

Needle Terminology

Length: Distance measured along the circumference of the needle, from the point to the swage. It is depicted on the packaging.

Shape: Different surgeries or surgical techniques demand different needle shapes.

Point: The sharpened part of the needle that first penetrates the tissue. The point consists of the area from the sharpened tip of the needle to the widest cross section of the needle body. The point determines ease of penetration and the initial size and shape of the hole.

Swage: Connection point of the suture and the needle. Micro teeth grip the suture material and prevent it from slipping out.

Wire Diameter: Gauge or thickness of the wire originally used to form the needle body.

Body: The part between the point and the swage of the needle.


Included Angle (Curvature): The fraction of the circle included in the curve of the needle. This is depicted on the packaging.

Needle Shapes And General Uses


Straight (AKA: KEITH): Skin Closure

3/8 Circle: Skin Closure, Ophthalmic


1/2 Circle: Skin Closure, GI, Cardiovascular, Ophthalmic

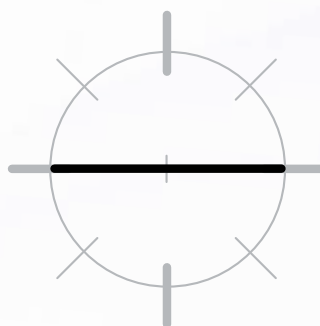

Straight

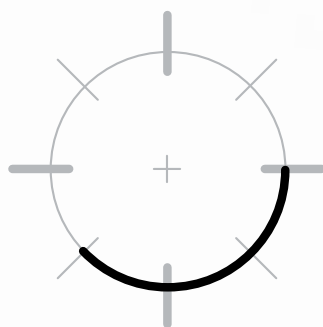



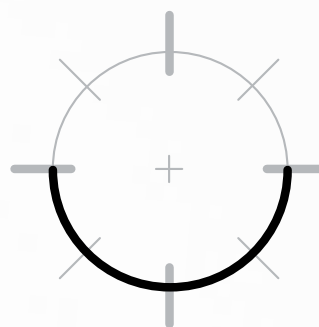





3/8 Circle



1/2 Circle

Needle Points And Cutting Edges

Needle Type	Point Geometry	Description	Commonly Used For
Conventional Cutting		Triangular point with cutting edge on the inner curvature	General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery
Precision Conventional Cutting		Triangular point with cutting edge on the inner curvature and precision-honed for added sharpness	General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery
Reverse Cutting		Triangular point with cutting edge on the outer curvature	Skin closure, ligamentous or fibrous tissues
Precision Reverse Cutting		Triangular point with cutting edge on the outer curvature and precision-honed for added sharpness	Primarily used in plastic surgery for delicate procedures and where a good cosmetic result is desired
Taper Point		Round shaft, straight or curved, taper point with no cutting edge	Soft tissue closure such as fascia, vascular, gastrointestinal, and most soft tissue below the skin
Taper Cutting		Four (4) cutting edges with a taper shaft	For penetration through tough tissue
Side Cutting Lancet		Spatula needle with the cutting edge on the inner, or outer curvature	Ophthalmic and microsurgery
Straight Cutting		Reverse cutting tip on a straight shaft	General skin closure

Needle Chart

Reverse Cutting Needles



LOOK™	Ethicon™	Syneture™
3/8 Circle		
C16	C-2	C-1 ————— 11mm
C17	C-2	C-1 ————— 12mm
•	C-3	C-12 ————— 14mm
C22	FS-3	C-12 ————— 16mm
C6	FS-2	C-13 ————— 18mm
C7	FS-1	C-14 ————— 24mm
C8	FS	C-15 ————— 26mm
C10	FSL	C-16 ————— 30mm
1/2 Circle		
C26	M-1	C-21 ————— 15mm
C19	•	• ————— 18mm
C20	J-1	C-22 ————— 19mm
C31	X-1	C-23 ————— 24mm
C12	CP-1	GS-11 ————— 27mm
C21	CP-2	GS-10 ————— 26mm
C13	CP-1	GS-11 ————— 37mm
C14	CP	GS-12 ————— 40mm

Precision Reverse Cutting Needles



LOOK™	Ethicon™	Syneture™
3/8 Circle		
C1	P-1	P-10 ————— 11mm
C3	P-3	P-13 ————— 13mm
PC34	PS-3	P-11 ————— 16mm
PC31	PS-2	P-12 ————— 18mm
PC33	PS-1	P-14 ————— 24mm

Precision Conventional Cutting Needles



LOOK™	Ethicon™	Syneture™
3/8 Circle		
CP1	PC-1	P-13 ————— 11mm

Needle Chart

Taper Point Needles



LOOK™	Ethicon™	Syneture™
1/2 Circle		
T22	TF	CV-22
T28	RB-1	CV-23
T14	SH-1	CV-24
T4	SH-1	CV-25
T	SH	V-20
T5	CT-2	GS-22
HR30	V-30	•
T9	MH	V-26
T8	CT-1	GS-21
T10	CT	GS-24

Straight Cutting Needles



LOOK™	Ethicon™	Syneture™
SC2	KS	SC-2

Straight Taper Point



LOOK™	Ethicon™	Syneture™
S	•	•

Taper Cutting Needles



LOOK™	Ethicon™	Syneture™
1/2 Circle		
Y5	V-5	KV-5
HRT31	V-34	KV-34
HRT36	V-34	KV-34

Side Cutting Lancet





LOOK™	Ethicon™	Syneture™
VSP8	S-24	SS-24

Plain Gut


Precision Conventional Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	1
CP1  11mm 3/8 Circle PC-1 P-13	Plain Gut	18"/45cm	1415B 						




Precision Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C3  13mm 3/8 Circle P-3 P-13	Plain Gut	18"/45cm		596B 686G SG5686	1241B 1644G SG5654G			

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C17  12mm 3/8 Circle C-2 C-1	Plain Gut	10"/25cm		521B 	531B 			
		18"/45cm		535B 686G CG172	536B 1644G SG721			
C6  18mm 3/8 Circle FS-2 C-13	Plain Gut	18"/45cm		500B  SG751	591B  SG752	592B  SG5630G		
		27"/70cm		550B H820G SG820G	551B H821H SG821G	554B H822H SG822G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).





• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Plain Gut

Reverse Cutting

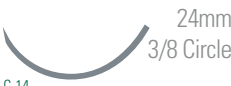

Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C7  FS-1 C-14	Plain Gut	27"/70cm				553B H810H SG810	555B H823H SG823	
C26  M-1 C-21	Plain Gut	18"/45cm			537B • •			
C31  X-1 C-23	Plain Gut	18"/45cm			599B • SG1633G	594B 612G SG612G		
C21  CP-2 GS-10	Plain Gut	18"/45cm				593B 612G SG612G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).




• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Chromic Gut

Precision Reverse Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
C1  11mm 3/8 Circle P-1 P-10	Chromic Gut	12"/30cm	1242B 796G G1780						
C3  13mm 3/8 Circle P-3 P-13	Chromic Gut	18"/45cm		1248B 687G SG5687G	595B 1654G SG5644G				
PC31  18mm 3/8 Circle PS-2 P-12	Chromic Gut	18"/45cm			1253B 1637G SG5637G				

Reverse Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
C16  11mm 3/8 Circle C-2 C-1	Chromic Gut	18"/45cm		543B 792G G1792K					
C17  12mm 3/8 Circle C-2 C-1	Chromic Gut	18"/45cm		563B 687G SG197	546B 1654G SG5644G				
C6  18mm 3/8 Circle FS-2 C-13	Chromic Gut	18"/45cm			558B 1637G SG675	525B 1638H SG676			
		27"/70cm		557B 634G SG634G	559B 635H SG635G	560B 636H SG636G			
C7  24mm 3/8 Circle FS-1 C-14	Chromic Gut	27"/70cm					561B 637H SG637G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).




• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Chromic Gut



Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C8  26mm 3/8 Circle FS C-15	Chromic Gut	27"/70cm						562B • •
C26  15mm 1/2 Circle M-1 C-21	Chromic Gut	18"/45cm		598B • •	544B 1653G •			
C20  20mm 1/2 Circle J-1 C-22	Chromic Gut	18"/45cm			547B 1643G •			

Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
T28  18mm 1/2 Circle RB-1 CV-23	Chromic Gut	14"/35cm		590B • UG202				
T4  22mm 1/2 Circle SH-1 CV-25	Chromic Gut	14"/35cm			589B G181H GG181			

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

PolySyn FA™ / Polyglycolic Acid


Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C6  18mm 3/8 Circle FS-2 C-13	Undyed Braided	27"/70cm				MC-423NC • •		
C7  24mm 3/8 Circle FS-1 C-14	Undyed Braided	27"/70cm				MC-442NC VC442H SL693		

Taper Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
HRT 36  36mm 1/2 Circle V-34 KV-34	Undyed Braided	36"/90cm					MC-917NC VR917 J9962T	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

PolySyn™/ Polyglycolic Acid

Precision Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C3 13mm 3/8 Circle P-3 P-13	Undyed Braided	18"/45cm	492B J492G SL5688G	493B J493G SL5687G	494B J494G SL5691G			
PC31 19mm 3/8 Circle PS-2 P-12	Undyed Braided	18"/45cm		495B J495G SL5636G	496B J496G SL5627G	497B J497G SL5628G		
PC33 24mm 3/8 Circle PS-1 P-14	Undyed Braided	18"/45cm				683B J683G SL5640G		

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C17 12mm 3/8 Circle C-2 C-1	Undyed Braided	10"/25cm		385B VC493G SL606	386B VC494G SL605			
C6 19mm 3/8 Circle FS-2 C-13	Violet Braided	27"/70cm			MC-397NC • •			
	Undyed Braided	27"/70cm		421B J421H SL690G	422B 422BC J422H SL691G	423B J423H SL822G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

PolySyn™ / Polyglycolic Acid


Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
C7 24mm 3/8 Circle  FS-1 C-14	Violet Braided	27"/70cm		GC-451NC • •				
	Undyed Braided	27"/70cm		441B J441H SL692G	GC-452NC 442B J442H SL693	443B J443H SL823		
C14 40mm 1/2 Circle  CP GS-12	Violet Braided	27"/70cm						480B J480H CL805

Reverse Cutting / Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C12/T8 27mm /36.6mm 1/2 Circle  CP-1 GS-11 CT-1 GS-21	Violet Braided	54"/135cm						GC-532NC J532H •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

PolySyn™/ Polyglycolic Acid









Taper Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
HRT 31  31mm 1/2 Circle V-34 KV-34	Violet Braided	30"/75cm							GC-936NC W9391 

Taper Point



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
T22  13mm 1/2 Circle TF CV-22	Undyed Braided	27"/70cm		GC-433NC  					
	Violet Braided	18"/45cm	GC-981NC  						
T28  18mm 1/2 Circle RB-1 CV-23	Violet Braided	30"/75cm			GC-304NC  				







Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

PolySyn™ / Polyglycolic Acid

Taper Point 

Needle	Material	Length	Length					
			5-0	4-0	3-0	2-0	0	1
T14  20mm 1/2 Circle SH-1 CV-24	Violet Braided	27"/70cm		310B/ 310BC J310H GL181	311B/ 311BC J311H GL182			
T5  26.50mm 1/2 Circle CT-2 GS-22	Violet Braided	27"/70cm			332B J332H CL882	333B J333H CL883		
T  26mm 1/2 Circle SH V-20	Violet Braided	27"/70cm			316B/ 316BC J316H GL122	317B/ 317BC J317H GL123		
HR30  30mm 1/2 Circle V-30	Violet Braided	30"/75cm			GC-322BC • •	GC-323NC • •		
T8  36.60mm 1/2 Circle	Violet Braided	27"/70cm				339B J339H CL930	340B J340H CL812	341B J341H CL813
		36"/90cm				GC-354NC • •	GC-518NC • •	
T10  40mm 1/2 Circle CT GS-24	Violet Braided	27"/70cm					352B J352H CL802	353B J353H CL803
		36"/90cm					GC-358NC • •	GC-359NC • •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).




• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Silk



Precision Reverse Cutting



Needle	Material	Length	Size					
			8-0	7-0	6-0	5-0	4-0	3-0
C1  11mm 3/8 Circle P-1 P-10	Black Braided	18"/45cm			1260B 639G SS1639G			
C3  13mm 3/8 Circle P-3 P-13	Black Braided	18"/45cm				1265B 640G SS5640G	792B 641G SS5641G	
PC31  18mm 3/8 Circle PS-2 P-12	Black Braided	18"/45cm					1270B 1677G SS5677G	

Reverse Cutting



Needle	Material	Length	Size					
			8-0	7-0	6-0	5-0	4-0	3-0
C16  11mm 3/8 Circle C-2 C-1	Black Braided	18"/45cm						795B • •
C17  12mm 3/8 Circle C-2 C-1	Black Braided	18"/45cm			752B • S1172	774B • S1173	787B 739G S1174	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Silk

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C22 16mm 3/8 Circle  FS-3 C-12	Black Braided	18"/45cm	750B • SS681					
C6 18mm 3/8 Circle  FS-2 C-13	Black Braided	18"/45cm		754B 682G SS682G	781B 683G SS683G	784B 1679H SS694		
		30"/75cm				768B 1679H SS745		
C7 24mm 3/8 Circle  FS-1 C-14	Black Braided	18"/45cm			776B 629G SS629	786B 684G SS684G		
		30"/75cm				775B • •		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).








• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Silk

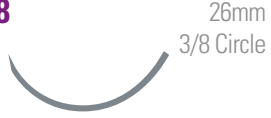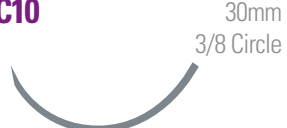Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C8  FS C-15	Black Braided	18"/45cm					770B 685G SS685G	
C10  FSL C-16	Black Braided	18"/45cm				765B 675H SS675	772B • •	
C26  M-1 C-21	Black Braided	18"/45cm		779B • •	777B M154T •	778B • •		
C19  	Black Braided	18"/45cm				783B • •		
C20  J-1 C-22	Black Braided	18"/45cm			780B 734G SS733			
C21  CP-2 GS-10	Black Braided	18"/45cm				782B • •		
C31  X-1 C-23	Black Braided	18"/45cm				785B 632G SS732		


Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.






LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Nylon

Precision Conventional Cutting

Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
CP1  11mm 3/8 Circle PC-1 P-13	Black Monofilament	18"/45cm	1456B 1856G SN1956		1454B 1854G SN1954				

Precision Reverse Cutting

Needle	Material	Length	Size						
			7-0	6-0	5-0	4-0	3-0	2-0	
C1  11mm 3/8 Circle P-1 P-10	Black Monofilament	18"/45cm		1279B 697G SN3697					
C3  13mm 3/8 Circle P-3 P-13	Black Monofilament	8"/20cm		1273B • •	1284B • •				
		18"/45cm		1282B 1698G SN5696G	1285B 698G SN5698G	1291B 699G SN5699G			
	Clear Monofilament	18"/45cm			1275B 690G SN5690G				
PC34  16mm 3/8 Circle PS-3 P-11	Black Monofilament	18"/45cm		1280B 1665G SN5767G	1283B 1668G SN5668G				
PC31  18mm 3/8 Circle PS-1 P-14	Black Monofilament	18"/45cm		1281B • SN5665G	1288B 1666G SN5666G	1290B 1667G SN5667G			
PC33  24mm 3/8 Circle PS-1 P-14	Black Monofilament	18"/45cm						1294B 1669G SN5669G	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).







• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Nylon

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C17 C-2 C-1  12mm 3/8 Circle	Black Monofilament	10"/25cm	911B • SN667G	913B • SN871	931B • SN772G			
	Clear Monofilament	10"/25cm		946B 668G SN5690G				
	Black Monofilament	18"/45cm	916B 667G SN667G	917B 668G SN871	921B 1954G SN872			
C22 FS-3 C-12  16mm 3/8 Circle	Black Monofilament	18"/45cm	915B 1665G SN660G	919B 1668G SN659G				
C6 FS-2 C-13  18mm 3/8 Circle	Black Monofilament	10"/25cm		945B 661H SN661G	944B 662G SN662G			
		18"/45cm		918B 661G SN661G	922B 662G SN662G	928B • SN653		
C7 FS-1 C-14  24mm 3/8 Circle	Black Monofilament	10"/25cm			947B 1629H SN658			
		18"/45cm			920B 1629H SN658	925B 663G SN663G		
C8 FS C-15  26mm 3/8 Circle	Black Monofilament	18"/45cm					927B 664G SN764G	
C10 FSL C-16  30mm 3/8 Circle	Clear Monofilament	27"/70cm				950B 1671H SN671		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).


• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Nylon


Straight Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
SC2 60mm Straight  KS SC-2	Black Monofilament	30"/75cm				923B 627H SN627		

Straight Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
S 61mm Straight 	Black Monofilament	30"/75cm					933B • •	


Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.



LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Polypropylene

Precision Conventional Cutting

Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
CP1  11mm 3/8 Circle PC-1 P-13	Blue Monofilament	8"/20cm	1057B • •	1058B • •				

Precision Reverse Cutting

Needle	Material	Length	Size					
			7-0	6-0	5-0	4-0	3-0	2-0
C1  11mm 3/8 Circle P-1 P-10	Blue Monofilament	18"/45cm		1079B 8697G SP5697G				
C3  13mm 3/8 Circle P-3 P-13	Blue Monofilament	8"/20cm		1077B • •	1078B • •			
		18"/45cm		1082B 8695G SP5695G	1085B 8698G SP5698G	1091B 8699G SP5699G		
PC31  18mm 3/8 Circle PS-2 P-12	Blue Monofilament	18"/45cm				1020B 8682G SP5682G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Polypropylene

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
C17  12mm 3/8 Circle C-2 C-1	Blue Monofilament	10"/25cm	1011B • •	1013B • •	1031B • •			
C6  18mm 3/8 Circle FS-2 C-13	Blue Monofilament	18"/45cm		1018B 8661G SP661G	1022B 8683G SP683G			
C7  24mm 3/8 Circle FS-1 C-14	Blue Monofilament	18"/45cm			1095B • •			

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).


• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Polyviolene™ / Polyester - Uncoated


Reverse Cutting



Needle	Material	Length	Size					
			9-0	8-0	7-0	6-0	5-0	4-0
C22  16mm 3/8 Circle FS-3 C-12	Green Braided Uncoated	18"/45cm					1176B X681G •	
C6  18mm 3/8 Circle FS-2 C-13	Green Braided Uncoated	18"/45cm						1180B X683G •
	White Braided Uncoated	18"/45cm						1205B R633G XTC692


Side Cutting Lancet



Needle	Material	Length	Size					
			9-0	8-0	7-0	6-0	5-0	4-0
VSP8  8mm 1/4 Circle S-24 SS-24	White Braided Uncoated	18"/45cm					FA-893NC W893 •	

Taper Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
Y5  16mm 1/2 Circle V-5 KV-5	Green Braided Uncoated	18"/45cm		797B X951H •				

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Additional Products


Silk Bulk Suture - Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
C7  24mm 3/8 Circle FS-1 C-14	Silk, Black Braided Strands	250 Pieces 18"/45cm			OF766 • •			
C8  26mm 3/8 Circle FS C-15	Silk, Black Braided Strands	250 Pieces 18"/45cm				OF770 • •		
C10  30mm 3/8 Circle FSL C-16	Silk, Black Braided Strands	250 Pieces 18"/45cm					OF772 • •	
C21  26mm 1/2 Circle CP-2 GS-10	Silk, Black Braided Strands	250 Pieces 18"/45cm			OF782 • •			

Silk Bulk Suture - Straight Cutting



NEEDLE	MATERIAL	LENGTH	SIZE					
			6-0	5-0	4-0	3-0	2-0	0
SC2  60mm Straight KS SC-2	Silk, Black Braided Strands	250 Pieces 30"/75cm				OF762 • •	OF767 • •	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Additional Products

Nylon Bulk Suture - Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
C6  18mm 3/8 Circle FS-2 C-13	Nylon, Black Monofilament Strands	250 Pieces 18"/45cm			OF928 • •			
C8  26mm 3/8 Circle FS C-15	Nylon, Black Monofilament Strands	250 Pieces 18"/45cm				OF927 • •		

LOOK™ Nylon Spools / Non-Sterile

Material	Length	Size					
		5-0	4-0	3-0	2-0	0	1
Non-Absorbable Nylon, Black Braided	100 yards	SP175 • •	SP176 • •				
Non-Absorbable Nylon, Black Monofilament	100 yards	SP185 • •	SP186 • •	SP187 • •	SP188 • •	SP189 • •	SP190 • •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Additional Products

LOOK™ Spools Polyviolene™/ Non-Sterile

Material	Length	Size					
		7-0	6-0	5-0	4-0	3-0	2-0
Non-Absorbable Polypropylene, Blue Monofilament	100 yards			SP135 • •	SP136 • •	SP137 • •	SP138 • •

LOOK™ Silk Spools / Non-Sterile

Material	Length	Size							
		6-0	5-0	4-0	3-0	2-0	0	1	2
Silk, Black Braided on Spools	25 yards	SP102 • •	SP103 • •	SP104 • •	SP105 • •	SP106 • •			
	50 yards								SP121 • •
	100 yards	SP114 • •	SP115 • •	SP116 • •	SP117 • •	SP118 • •	SP119 • •	SP120 • •	
Silk, White Braided on Spools	2000 yards				SP1230 • •				

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

•• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Additional Products

LOOK™ Polyviolene™ Spools / Non-Sterile

Material	Length	Size					
		5-0	4-0	3-0	2-0	0	1
Non-Absorbable Polyester, Green Braided, Coated	250 yards		SP252/ SP252R ● ●				
	3281 yards		SP2520 ● ●				

LOOK™ Polyviolene™ Spools / Non-Sterile

Material	Length	Size					
		0	1	2	3	4	5
Non-Absorbable Polyester, Green Braided, Coated	100 yards	SP255 ● ●		SP257 ● ●			
	300 yards						SP2770 ● ●
Non-Absorbable Polyester, White Braided, Coated	100 yards			SP267 ● ●			

LOOK™ Polyviolene™ Spools, Uncoated / Non-Sterile

Material	Length	Size						
		5-0	4-0	3-0	2-0	0	1	2
Non-Absorbable Polyester, Green Braided, Uncoated	100 yards	SP201 ● ●		SP203 ● ●	SP204 ● ●		SP206 ● ●	SP207 ● ●
Non-Absorbable Polyester, White Braided, Uncoated	100 yards	SP211 ● ●	SP210 ● ●					

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

● indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

Product Index

Code	Size	Material	Page Number
0F762	3-0	Silk, Black Braided Strands	27
0F766	3-0	Silk, Black Braided Strands	27
0F767	2-0	Silk, Black Braided Strands	27
0F770	2-0	Silk, Black Braided Strands	27
0F772	0	Silk, Black Braided Strands	27
0F782	3-0	Silk, Black Braided Strands	27
0F927	2-0	Nylon, Black Monofilament Strands	28
0F928	3-0	Nylon, Black Monofilament Strands	28
1011B	6-0	Polypropylene, Blue Monofilament	25
1013B	5-0	Polypropylene, Blue Monofilament	25
1018B	5-0	Polypropylene, Blue Monofilament	25
1020B	4-0	Polypropylene, Blue Monofilament	24
1022B	4-0	Polypropylene, Blue Monofilament	25
1031B	4-0	Polypropylene, Blue Monofilament	25
1058B	5-0	Polypropylene, Blue Monofilament	24
1075B	6-0	Polypropylene, Blue Monofilament	24
1077B	6-0	Polypropylene, Blue Monofilament	24
1078B	5-0	Polypropylene, Blue Monofilament	24
1079B	6-0	Polypropylene, Blue Monofilament	24
1082B	6-0	Polypropylene, Blue Monofilament	24
1085B	5-0	Polypropylene, Blue Monofilament	24
1091B	4-0	Polypropylene, Blue Monofilament	24
1095B	4-0	Polypropylene, Blue Monofilament	25
1176B	5-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
1180B	4-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
1205B	4-0	Polyviolene™ / Polyester, White Braided Uncoated	26
1241B	4-0	Plain Gut	9
1242B	6-0	Chromic Gut	11
1248B	5-0	Chromic Gut	11
1253B	4-0	Chromic Gut	11
1260B	6-0	Silk, Black Braided	18
1265B	5-0	Silk, Black Braided	18
1270B	4-0	Silk, Black Braided	18
1273B	6-0	Nylon, Black Monofilament	21
1275B	5-0	Nylon, Clear Monofilament	21
1276B	7-0	Nylon, Black Monofilament	21
1279B	6-0	Nylon, Black Monofilament	21
1280B	6-0	Nylon, Black Monofilament	21
1281B	6-0	Nylon, Black Monofilament	21
1282B	6-0	Nylon, Black Monofilament	21

Product Index

Code	Size	Material	Page Number
1283B	5-0	Nylon, Black Monofilament	21
1284B	5-0	Nylon, Black Monofilament	21
1285B	5-0	Nylon, Black Monofilament	21
1288B	5-0	Nylon, Black Monofilament	21
1290B	4-0	Nylon, Black Monofilament	21
1291B	4-0	Nylon, Black Monofilament	21
1294B	3-0	Nylon, Black Monofilament	21
1415B	6-0	Plain Gut	9
1454B	4-0	Nylon, Black Monofilament	21
1456B	6-0	Nylon, Black Monofilament	21
310B	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
310BC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
311B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
311BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
316B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
316BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
317B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
317BC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
332B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
333B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
339B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
340B	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
341B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
352B	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
353B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
385B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
386B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
421B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
422B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
422BC	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
423B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
441B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
442B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
443B	2-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
480B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	15
492B	6-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
493B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
494B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
495B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
496B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14

Product Index

Code	Size	Material	Page Number
497B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
500B	5-0	Plain Gut	9
521B	5-0	Plain Gut	9
525B	3-0	Chromic Gut	11
531B	4-0	Plain Gut	9
535B	5-0	Plain Gut	9
536B	4-0	Plain Gut	9
537B	4-0	Plain Gut	10
543B	5-0	Chromic Gut	11
544B	4-0	Chromic Gut	12
547B	4-0	Chromic Gut	12
550B	5-0	Plain Gut	9
551B	4-0	Plain Gut	9
553B	3-0	Plain Gut	10
554B	3-0	Plain Gut	9
555B	2-0	Plain Gut	10
557B	5-0	Chromic Gut	11
558B	4-0	Chromic Gut	11
559B	4-0	Chromic Gut	11
560B	3-0	Chromic Gut	11
561B	2-0	Chromic Gut	11
562B	0	Chromic Gut	12
563B	5-0	Chromic Gut	11
564B	4-0	Chromic Gut	11
589B	4-0	Chromic Gut	12
590B	5-0	Chromic Gut	12
591B	4-0	Plain Gut	9
592B	3-0	Plain Gut	9
593B	3-0	Plain Gut	10
594B	3-0	Plain Gut	10
595B	4-0	Chromic Gut	11
596B	5-0	Plain Gut	9
598B	5-0	Chromic Gut	12
599B	4-0	Plain Gut	10
683B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
750B	6-0	Silk, Black Braided	19
752B	6-0	Silk, Black Braided	18
754B	5-0	Silk, Black Braided	19
765B	3-0	Silk, Black Braided	20
768B	3-0	Silk, Black Braided	19

Product Index

Code	Size	Material	Page Number
770B	2-0	Silk, Black Braided	20
772B	0	Silk, Black Braided	20
774B	5-0	Silk, Black Braided	18
775B	3-0	Silk, Black Braided	19
776B	4-0	Silk, Black Braided	19
777B	4-0	Silk, Black Braided	20
778B	3-0	Silk, Black Braided	20
779B	5-0	Silk, Black Braided	20
780B	4-0	Silk, Black Braided	20
781B	4-0	Silk, Black Braided	19
782B	3-0	Silk, Black Braided	20
783B	3-0	Silk, Black Braided	20
784B	3-0	Silk, Black Braided	19
785B	3-0	Silk, Black Braided	20
786B	3-0	Silk, Black Braided	19
787B	4-0	Silk, Black Braided	18
792B	4-0	Silk, Black Braided	18
795B	3-0	Silk, Black Braided	18
797B	5-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
911B	6-0	Nylon, Black Monofilament	22
913B	5-0	Nylon, Black Monofilament	22
915B	6-0	Nylon, Black Monofilament	22
916B	6-0	Nylon, Black Monofilament	22
917B	5-0	Nylon, Black Monofilament	22
918B	5-0	Nylon, Black Monofilament	22
919B	5-0	Nylon, Black Monofilament	22
920B	4-0	Nylon, Black Monofilament	22
921B	4-0	Nylon, Black Monofilament	22
922B	4-0	Nylon, Black Monofilament	22
923B	3-0	Nylon, Black Monofilament	23
925B	3-0	Nylon, Black Monofilament	22
927B	2-0	Nylon, Black Monofilament	22
928B	3-0	Nylon, Black Monofilament	22
931B	4-0	Nylon, Black Monofilament	22
933B	2-0	Nylon, Black Monofilament	23
944B	4-0	Nylon, Black Monofilament	22
945B	5-0	Nylon, Black Monofilament	22
946B	5-0	Nylon, Clear Monofilament	22
947B	4-0	Nylon, Black Monofilament	22
950B	3-0	Nylon, Clear Monofilament	22

Product Index

Code	Size	Material	Page Number
FA-893NC	5-0	Polyviolene™ / Polyester, White Braided Uncoated	26
GC-304NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
GC-322BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-323NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-354NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-358NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-359NC	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-433NC	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	16
GC-451NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	15
GC-452NC	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
GC-518NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-532NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	15
GC-936NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
GC-981NC	6-0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
MC-397NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	14
MC-423NC	3-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
MC-442NC	3-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
MC-917NC	2-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
SP102	6-0	Silk, Black Braided Spool, 25 yards	29
SP103	5-0	Silk, Black Braided Spool, 25 yards	29
SP103	4-0	Silk, Black Braided Spool, 25 yards	29
SP105	3-0	Silk, Black Braided Spool, 25 yards	29
SP106	2-0	Silk, Black Braided Spool, 25 yards	29
SP114	6-0	Silk, Black Braided Spool, 100 yards	29
SP115	5-0	Silk, Black Braided Spool, 100 yards	29
SP116	4-0	Silk, Black Braided Spool, 100 yards	29
SP117	3-0	Silk, Black Braided Spool, 100 yards	29
SP118	2-0	Silk, Black Braided Spool, 100 yards	29
SP119	0	Silk, Black Braided Spool, 100 yards	29
SP120	1	Silk, Black Braided Spool, 100 yards	29
SP121	2	Silk, Black Braided Spool, 50 yards	29
SP1230	3-0	Silk, Black Braided Spool, 2000 yards	29
SP135	5-0	Polypropylene, Blue Monofilament, 100 yards	29
SP136	4-0	Polypropylene, Blue Monofilament, 100 yards	29
SP137	3-0	Polypropylene, Blue Monofilament, 100 yards	29
SP138	2-0	Polypropylene, Blue Monofilament, 100 yards	29
SP175	5-0	Nylon, Black Braided Spool, 100 yards	28
SP176	4-0	Nylon, Black Braided Spool, 100 yards	28
SP185	5-0	Nylon, Black Monofilament Spool, 100 yards	28
SP186	4-0	Nylon, Black Monofilament Spool, 100 yards	28

Product Index

Code	Size	Material	Page Number
SP187	3-0	Nylon, Black Monofilament Spool, 100 yards	28
SP188	2-0	Nylon, Black Monofilament Spool, 100 yards	28
SP189	0	Nylon, Black Monofilament Spool, 100 yards	28
SP190	1	Nylon, Black Monofilament Spool, 100 yards	28
SP201	5-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP203	3-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP204	2-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP206	1	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP207	2	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP210	4-0	Polyviolene™ / Polyester, White Braided, Uncoated, Spool, 100 yards	30
SP211	5-0	Polyviolene™ / Polyester, White Braided, Uncoated, Spool, 100 yards	30
SP252	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 250 yards	30
SP2520	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 3281 yards	30
SP252R	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 250 yards	30
SP2555	0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 100 yards	30
SP257	2	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 100 yards	30
SP267	2	Polyviolene™ / Polyester, White Braided, Coated, Spool, 100 yards	30
SP2770	5	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 300 yards	30

Visit SurgicalSpecialties.com And Click "Find A Product" For SKUs And Cross-Reference Information.

Surgical Specialties Corporation, 50 Braintree Hill Office Park, Suite 101, Braintree, MA 02184, U.S.A.
Phone: 610.404.3520 • 877.991.1110 (USA and Canada only) • Fax: 610.404.4010
All trademarks are property of their respective owners. ©2015 Surgical Specialties Corporation.
All Rights Reserved. SSC-027 R0 5/15

**SURGICAL
SPECIALTIES**
CORPORATION™