


3M™ Steri-Drape™ Shoulder Split Sheet with Pouch – 9196

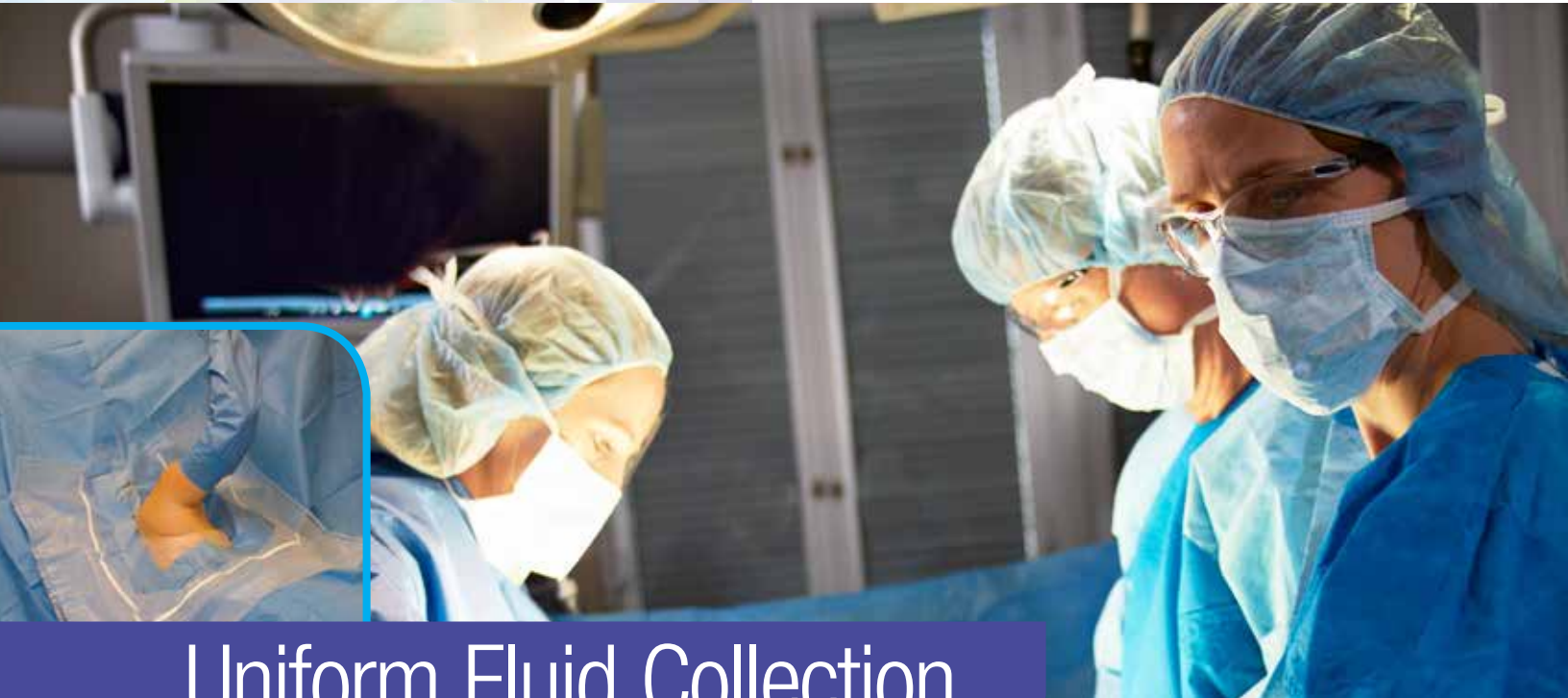

Uniform Fluid Collection

Reduced Exposure Risk

3M Steri-Drape Shoulder Split Sheet has a wide, rectangular fluid collection pouch that collects fluids uniformly, which can be positioned above or below the shoulder. Paired with 3M Fabric, it helps reduce the risk of exposure to blood and body fluids for patients and healthcare workers.

- Designed for one-piece draping in lateral or beach chair position
- Reliable 2" adhesive strip allows the drape to adhere better throughout the procedure
- Adhesive liner is separated in four pieces for added control during application
- Exit ports on both sides of pouch help deliver better fluid control
- Malleable metal wire keeps the pouch open for optimal fluid collection
- The fabric absorbs spills and splatters, making multiple layers unnecessary and prevents strike-through


Lead
The Way

3M

3M™ Steri-Drape™ Shoulder Split Sheet with Pouch – 9196

3M Fabric

- Meets Flammability Class 1 (CFR Title 16 Part 1610)
- Meets AAMI Standard PB70:2012 *Liquid barrier performance and classification of protective apparel and drapes intended for use in health care facilities, Selection of surgical gowns and drapes in health care facilities*
- Provides AAMI level 4 barrier performance throughout entire drape
- Low linting (INDA IST 160 1-83)
- Meets AORN 2014 Recommended Practices for materials, selection, and use of surgical drapes that represent an optimal level of practice

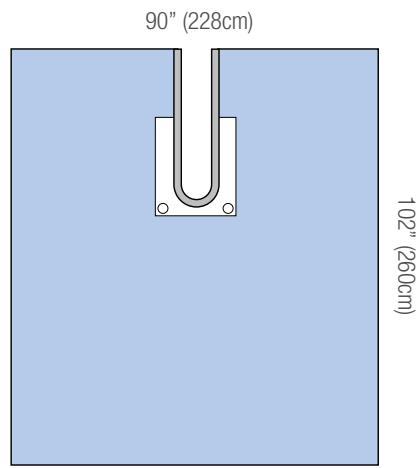

Ordering Information

3M™ Steri-Drape™ Shoulder Split Sheet with Pouch – 9196

Cat No.	Product Description	Packaging	Minimum Order	Pricing Unit
9196	90" x 102" (228cm x 260cm)	1 Each 10 Each/Box 2 Box/Case	Case	Box

For more information about 3M™ Surgical Drapes or other 3M Health Care products:

- » Contact your local 3M representative
- » Call the 3M Health Care Customer Helpline at 800-228-3957
- » Visit www.3M.com/infectionprevention



Infection Prevention Division 3M Health Care

2510 Conway Avenue
St. Paul, MN 55144-1000
U.S.A.
1-800-228-3957
www.3M.com/infectionprevention

3M and Steri-Drape are trademarks of 3M.
©3M 2014. All rights reserved.
70-2009-4653-4