


# CFUR301 & CFUR401 *Urology Tables*

- 53" metal-free imaging area
- 3 table extensions included
- 500 lb patient capacity
- 3 or 4 motion models



ISO 13485:2003  
FM 584901

made in the USA  
with US & imported parts

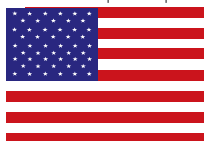

OAKWORKS<sup>®</sup>  
**MEDICAL**  
*Position for Success*

# CFUR301 & CFUR401 Urology Tables

## 3 Radiolucent Extensions

### 33" WIDE SUPPORT SURFACE

The two lateral extensions give you the support area you need while maximizing access to the 23" W x 47" L Carbon Fiber Top. Imaging Area 19" W x 24" L (48cm x 61cm)

*Radiolucent patient transfer extension creates 53" long metal-free imaging area between the T-rails.*



### (2) 5" X 24" RADIOLUCENT LATERAL EXTENSIONS #66158



Removable lateral extensions slide and lock on T-rails.

### (1) 23" WIDE X 28" LONG RADIOLUCENT PATIENT TRANSFER EXTENSION #66156



Lightweight design can be used for general imaging when a full length diving board table design is required.

## CFUR301 #75214-T01      CFUR401 #75215-T01

STANDARD		OPTIONS <i>must be ordered with table.</i>		#
Carbon Fiber Top 53" (135cm) Metal free imaging space with Patient Transfer Extension	23" (58cm) wide x 47" (119cm) long	5" (13cm) casters	upgrade from the standard 4" (10cm)	65427
<b>CFUR301 MOVEMENTS</b>		Foot Control (CFUR301)	multi-functional	75076
Motion 1	Electric Height Range 25.5" - 43.5" (65-110cm)	Foot Control (CFUR401)	multi-functional	75078
Motion 2	Electric Trendelenburg tilt +15°/-11°	Non-USA electric system	please contact us for order details	73134
Motion 3	Electric Longitudinal travel 10" (25cm)	<b>ACCESSORIES <i>order any time.</i></b>		#
<b>CFUR401 MOVEMENTS</b>		Protective Drape	Clear PVC - attaches to T-Rail	69390
Motion 1	Electric Height Range 26.5" - 44.5" (67-114cm)	Radiolucent Lateral Extension	5" x 24" (13cm x 61cm)	66158
Motion 2	Electric Trendelenburg tilt +15°/-11°	Radiolucent Lateral Extension	5" x 14" (13cm x 36cm)	69794
Motion 3	Electric Longitudinal travel 10" (25cm)			
Motion 4	Electric Lateral travel ±4" (10cm)			
Controls	Multi-function hand control			
Weight Capacity	500 lb (227kg) Each Extension: 200 lb (91kg)			
Product Weight	CFUR301 568 lb (258kg) CFUR401 590 lb (268kg)			
Radiolucent Lateral Extension	(2) 5" x 24" (12cm x 61cm)			
Radiolucent Patient Transfer Extension	(1) 23" x 28" (58cm x 71cm)			
Included Features	1" (3cm) Radiolucent Comfort Foam™ table top pad, 4" (10cm) locking casters, Digital angle gauge, patient safety straps (2), Industry standard full length T-rails on each side, Urology bag hoop			
Warranty	3 years parts, 2 years labor			
Safety Listings	FDA listed, CE marked, cETLus approved, CB Certificate			

### UROLOGY BAG HOOP

Table comes with hoop and clamps to t-rails. (bags not included)



INCLUDED

### 5" CASTER #65427

Recommended for mobile service providers



OPTIONAL

# UROLOGY OPTIONS & ACCESSORIES

## PROTECTIVE DRAPE

#69390

This removable clear vinyl drape acts as a splash shield. It attaches to the t-rail and facilitates movement of C arms and lithotripters without concerns about bumping into the table mechanisms.



ACCESSORY

## ADJUSTABLE IV POLE WITH INTEGRATED CLAMP

#72871



ACCESSORY

## LIFT ASSIST STIRRUP SYSTEM

*BARIATRIC RATED* #72868



ACCESSORY

## STAINLESS STEEL

CLARK SOCKET #72870



ACCESSORY

## ADJUSTABLE KNEE

CRUTCH SUPPORT #72872



ACCESSORY

## LIFT ASSIST STIRRUP


SYSTEM #72866



ACCESSORY

## TOTAL KNEE

STABILIZER #72873

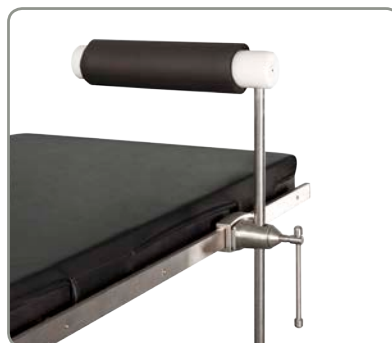

ACCESSORY

## BASIC STIRRUP

LEG HOLDERS #72869



ACCESSORY

# OAKWORKS<sup>®</sup> **MEDICAL** *Position for Success*

## CFUR301 & CFUR401 UROLOGY TABLES



### INCLUDED FEATURES:

**(3) Radiolucent Extensions**

**(2) Patient Safety Straps**  
84" L with 24" of Velcro<sup>®</sup> at each end



**1" Comfort Foam<sup>™</sup>  
Table Top Pad**



**4" Locking  
Casters**

**Industry Standard Full  
Length T-rails on each side**

The Integrated T-rails are US standard size of .375" (9.5mm) thick by 1.25" (28.6mm) wide.

**Digital Angle Gauge**

**Urology Bag Hoop**

Model: Urocatch<sup>®</sup> with intergrak clamps  
(Bag not included)



ISO 13485:2003  
FM 584901

Oakworks Medical is ISO 13485:2003 certified. The ISO standard is intended to ensure that medical device manufacturers have the required quality management systems in place to safely design, manufacture and distribute medical devices. This registration certifies that Oakworks Medical's quality system conforms to these requirements in our state-of-the-art 90,000 square foot factory in Pennsylvania.

CONTACT:

*View our full line of  
Medical products*



Scan this code with  
your smart phone

#72703



717.235.6807 [www.oakworksmed.com](http://www.oakworksmed.com) email: [customerservice@oakworks.com](mailto:customerservice@oakworks.com)

65734 REVISED 08.25.15