# Premier

#### **IMPORTANT NOTE:**

Each test can be used only ONE TIME. Do not attempt to use the test more than once.

#### **Test Procedure**

### COVID-19 IgG And IgM Rapid Test Cassette Quick Reference Guide
































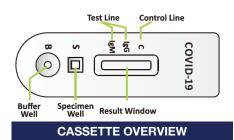



www.premierbiotech.com

### **POSITIVE RESULTS**



**NEGATIVE RESULT** 

IgG and IgM Antibodies Positive **IgG Antibodies Positive** Three lines appear(C, IgG and IgM) Two lines appear(C and IgG)

IgM Antibodies Positive Two lines appear(C and IgM)

## 

## 


**INCONCLUSIVE RESULTS** 

Review the procedure and repeat the test

Control line fails to appear.

with a new test cassette.



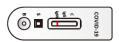



To learn more and to view our guick COVID-19 procedure video, visit our support page at: https://premierbiotech.com/innovation/covid-19/

Only one colored line appears

in the control region (C).