




OAKWORKS

Professional Massage Equipment

The #1 choice of massage therapists for over 35 years!

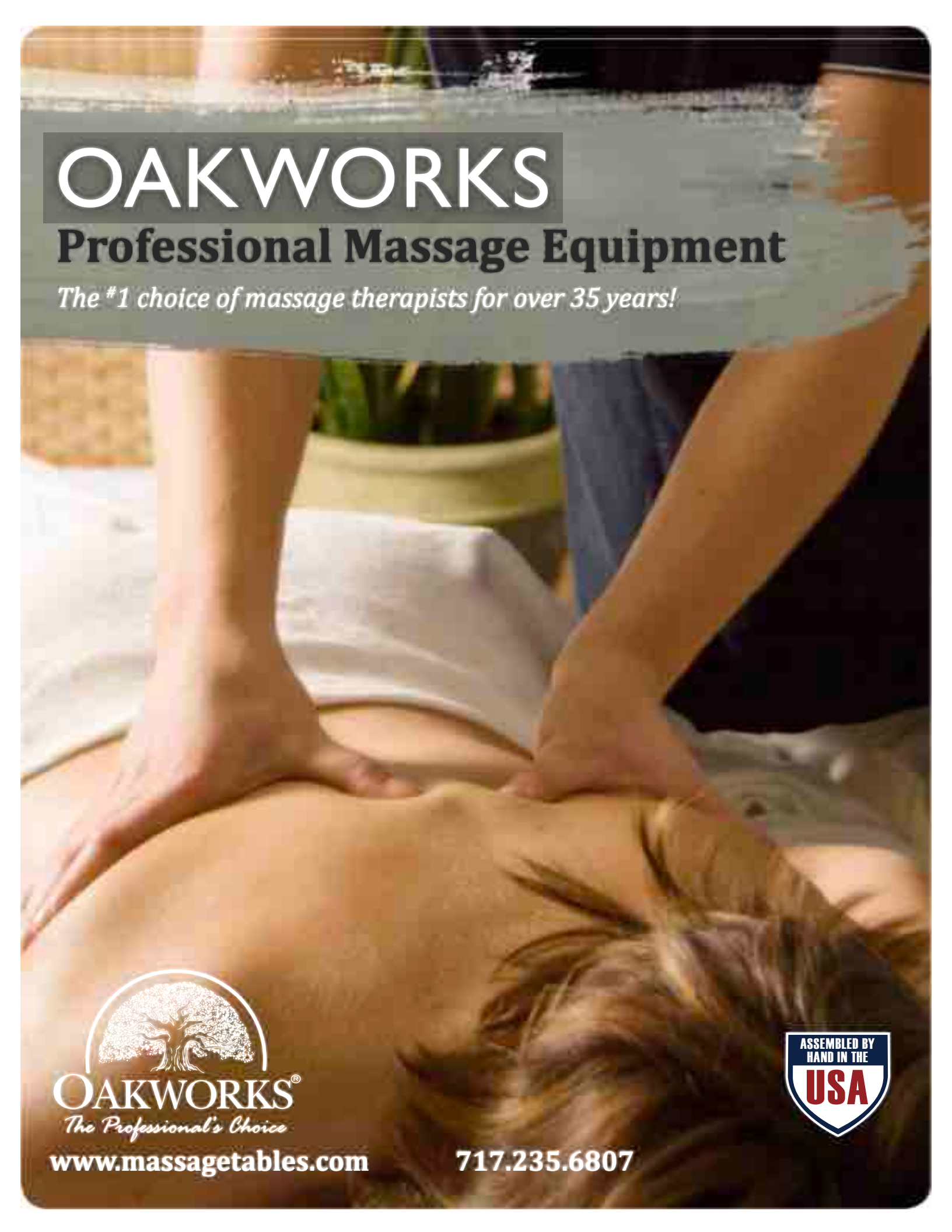

www.massagetables.com

717.235.6807



TABLE OF CONTENTS

PORTABLE TABLES & CHAIRS

Nova™ Table.....	4
WellSpring™ Table.....	5
Alliance™ Table.....	6
Aurora™ Table.....	7
Advanta™ Table	8
One Table.....	9
Portable Table Packages	10-11
Portal Pro® Chair.....	12
Portal Light Chair	13
Storable Mat	14

STATIONARY TABLES

ProLuxe™ Series.....	15
ProLuxe™ Convertible	16-17
New! PT400M.....	18-19
Clinician Series	20-21

ACCESSORIES

Boiance™ Face Rest.....	22
New! Boiance™ Float.....	22
Bolsters & Face Rests	23
Covers & Cases.....	24
Side Lying Positioning System.....	25
Stools & Table Cart.....	26

FABRIC & PADDING

TerraTouch™	27
Padding.....	27



*"Professional Massage Equipment
The Essentials For Success"*



Watch the video!

visit www.youtube.com/oakworksinc for more videos

ABOUT US

MADE IN THE USA



Jeff & Linda Riach

Since 1977, when my wife Linda and I made all of our tables by hand in the commune where we started Oakworks®, we have always understood the importance of US manufacturing. Back then, Oakworks® allowed us to work close to nature, with our little solar heated shop and organic gardens outside our door. We spent seven years living off the land and on the land, had two children, and developed a view of life and values that were different than most. We came to understand that we are all connected to the land and to each other, we have a responsibility to nurture all. This sense of connection and responsibility is reflected in our commitment to US manufacturing, the environment, our customers, and our employees.

AWARD WINNING CUSTOMER SERVICE

Our commitment to the environment and our employees is carried through to our commitment to our customers. On the rare occasion of a problem, our award winning Customer Service team is available to assist. We are committed to our customers' satisfaction, and want to be sure you are pleased with every purchase. Our reputation for Customer Service is the best in the industries we serve.



www.massagetables.com

923 East Wellspring Rd. New Freedom, PA 17349

DESIGN PHILOSOPHY

Our design philosophy is really quite simple; build beautiful products that support therapist's ergonomic health and client's comfort while ensuring good value through high quality construction. Beauty is in the eye of the beholder, and within this catalog we are sure you will see something that catches your eye. Ergonomics are critical to insure that your therapists deliver their best work, day after day, all day long. High quality construction is both environmental and customer imperative. We strive to deliver the highest quality products built to withstand the demands of any treatment room. We know that products that are built to last don't end up in landfills and we build products that are designed to be repaired if necessary, not thrown away.

RESPONSIBLE FOREST MANAGEMENT

Since our inception in 1977, we have maintained a strong commitment to both the environment and our employees. We make products that are built to last while having the least amount of impact on the environment as possible. We are Global Cooling Partners with Trees for the Future and have financed the planting of over 375,000 trees, completely offsetting our carbon footprint in recent years. We are a Forest Stewardship Council™ Certified plant, following their high standards that ensure we use woods harvested in responsible ways. Our Pennsylvania manufacturing facility meets or exceeds all EPA and OSHA standards for safe working conditions and we have developed many unique ergonomic work stations to insure our employees' health.



The mark of
responsible forestry

NOVA™

professional line

Constructed from the finest materials in the industry, the Nova™ is the most customizable and reliable table on the market. The Nova™ boasts an exceptional variety of options, widths and height ranges to easily accommodate any size therapist and/or bodywork modality. A table built to last a lifetime, the Nova™ has been the preferred choice among massage professionals for over 25 years.



Shown in Espresso with optional QuickLock™ face rest platform & Aero•Cel™ face rest crescent



1. **PATENTED SLIDELOCK™**
Designed to offer even greater strength than our cable lock system, features the strongest cables on the market.
2. **PATENTED INTEGRAHINGE™**
OAKWORKS® exclusive design integrates the attachment of the bracing system to offer incredible strength and lighter weight that no piano hinge can match.
3. **ENGRAVED WOODEN KNOBS**
Low profile design with built in grip band is easy to use and won't get in your way like plastic knobs.
4. **UNILOCK™**
Super quick, easy to use system folds out of the way for shiatsu and won't rip your sheets or mark your floor.
5. **ULTIMATE ACCESS END PANEL & SHIATSU CABLES**
Our unique system allows you to get your legs under the table at the height you work, offering the best access on the market. The shiatsu cables allow the table to be laid flat on the floor.
6. **AERO•CEL™ PADDING UPGRADE**
New technology in padding options now gives you 4" of more resilient support with a softness that conforms to the body's curves.

specifications

MAX. OPERATING WEIGHT	550 lbs.
WEIGHT	29-36 lbs
HEIGHT RANGES	18"-26", 22"-30", 24"-34", 26"-36"
WIDTHS	27", 29", 31", 33"
LENGTH	73"
TOPS	Rectangular, Round Corner (not available on 33" width)
PADDING (CFC FREE)	Plush, Aero•Cel™, Semi-Firm, Firm
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Ultimate Access End Panel, Shiatsu Cables
ACCESSORIES	Wooden Portable Arm Rest *, QuickLock™ Face Rest Platform, Face Hole w/Filler, and our full line of accessories
Bold - Standard • Italic - Additional Charge * 24"-34" and 26"-36" height ranges only	

UniLock™: U.S. Patent # 7,472,656

WELLSPRING™

professional line



Our lightest, full size table features our patented complementary support system and the best access on the market. This stylish lightweight table is supported by a sturdy ultimate access aluminum understructure ensuring superior leg room for the therapist on all sides of the table. This is a perfect table for anyone who travels, treats a diverse client base, or works seated much of the time.



Shown in Sapphire with optional QuickLock™ face rest platform & Aero-Cel™ face rest crescent



specifications

MAX. OPERATING WEIGHT	550 lbs.
WEIGHT	25-30 lbs
HEIGHT RANGES	17"-26", 22"-31", 24"-34"
WIDTHS	27", 29", 31"
LENGTH	73"
TOPS	Rectangular , <i>Rounded Corners</i>
PADDING (CFC FREE)	Semi-Firm , <i>Aero-Cel, Plush, Firm</i>
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Ultimate Access Ends
ACCESSORIES	<i>QuickLock™ Face Rest Platform, Arm Hammock, Field Feet, and our full line of accessories</i>
Bold - Standard • <i>Italic - Additional Charge</i>	

- PATENTED INTEGRAHINGE™**
OAKWORKS® exclusive design integrates the attachment of the bracing system to offer incredible strength and lighter weight that no piano hinge can match.
- UNILOCK™**
Super quick, easy to use closure system. Won't rip your sheets.
- COMPLEMENTARY SUSPENSION SYSTEM**
Designed to be the strongest system on the market while ensuring the best access to your client by allowing you to get your legs under the sides and the ends of the table, at any height!
- EXCEPTIONAL FRAME STRENGTH**
Frames are built on a computerized machine for accuracy and strength, utilizing box joints that are glued and screwed together.
- TELESCOPING LEGS**
Easy push button adjustment for quick height changes.
- AERO-CEL™ PADDING UPGRADE**
New technology in padding options now gives you 4" of more resilient support with a softness that conforms to the body's curves.

CS System™ (Complementary Suspension) : U.S. Patent #6,192,809
UniLock™: U.S. Patent # 7,472,656

ALLIANCE™ *professional line*

wood

The Alliance™ is the most versatile wood portable table on the market today. This table is perfect for any massage therapist with a need to perform any of today's most popular treatments such as esthetics, spa treatments and massage.

Shown in Ruby with optional QuickLock™ face rest platform & Aero•Cel™ face rest crescent



aluminum

The Alliance™ Aluminum is a lightweight table that offers exceptional strength, unbeatable client access, and positioning accuracy. This table is perfect for any of today's most popular treatments in the spa, medi-spa, and massage fields.



Shown in Orchid with optional QuickLock™ face rest platform & Aero•Cel™ face rest crescent

1. PATENTED INTEGRAL HINGE™

Oakworks® exclusive design integrates the attachment of the bracing system to offer incredible strength and lighter weight that no piano hinge can match.

2. UNILOCK™

Super quick, easy to use closure system. Won't rip your sheets.

3. ULTIMATE ACCESS END PANEL & SHIATSU CABLES (WOOD ONLY)

Our unique system allows you to get your legs under the table at the height you work, offering the best access on the market. The shiatsu cables allow the table to be laid flat on the floor.

4. EXCEPTIONAL FRAME STRENGTH

Frames are built on a computerized machine for accuracy and strength, utilizing box joints that are glued and screwed together.



specifications

MAX. OPERATING WEIGHT	550 lbs.
WEIGHT	32-38 lbs (wood) 30-36 lbs (aluminum)
HEIGHT RANGE	24"-34"
WIDTHS	27", 30"
LENGTH	73"
TOP	Back Rest
PADDING (CFC FREE)	Semi-Firm, Plush, Firm
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Ultimate Access End Panel, Shiatsu Cables - (wood only)
ACCESSORIES	QuickLock™ Face Rest Platform, Arm Hammock, and our full line of accessories
Bold - Standard • <i>Italic - Additional Charge</i>	

CS System™ (Complementary Suspension) : U.S. Patent #6,192,809
UniLock™: U.S. Patent # 7,472,656

AURORA™

professional line

A lightweight wood table specifically constructed for the massage professional whose practice requires constant travel. Designed to withstand the rigors of everyday use, the Aurora™ offers ultimate client comfort, and the strength to last an entire career. Truly a flawless investment for any massage therapist.



CE

specifications

- PATENTED INTEGRAHINGE™**
OAKWORKS® exclusive design integrates the attachment of the bracing system to offer incredible strength and lighter weight that no piano hinge can match.
- ENGRAVED WOODEN KNOBS**
Low profile design with built in grip band is easy to use and won't get in your way like plastic knobs.
- UNILOCK™**
Super quick, easy to use closure system folds that out of the way for shiatsu and won't rip your sheets or mark your floor.
- ULTIMATE ACCESS END PANEL & SHIATSU CABLES**
Our unique system allows you to get your legs under the table at the height you work, offering the best access on the market. The shiatsu cables allow the table to be laid flat on the floor.
- CABLE LOCK SYSTEM**
Our system uses stronger cables than anyone else, designed not to stretch or fail even under the heaviest loads.

MAX. OPERATING WEIGHT	550 lbs.
WEIGHT	25-31 lbs
HEIGHT RANGE	24"-34"
WIDTHS	28", 30"
LENGTH	73"
TOP	Rectangular
PADDING (CFC FREE)	Plush, Sem-Firm
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Access End Panel, Shiatsu Cables
ACCESSORIES	QuickLock™ Face Rest Platform, Wooden Portable Arm Rest, and our full line of accessories
Bold - Standard • Italic - Additional Charge	

UniLock™: U.S. Patent # 7,472,656

ADVANTA™ *professional line*




Shown in Ruby with optional QuickLock™ face rest platform & Aero•Cel™ face rest crescent

Optional Recess Fillers:
Attach as needed for additional width.
Shown in Opal.




The Advanta™ is the ultimate lightweight table for those therapists looking to maximize proper ergonomics. This unique table incorporates a special contour top that allows for superior therapist access to any size client as well as providing unmatched leg room underneath the table.

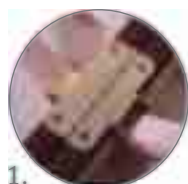

1.




2.

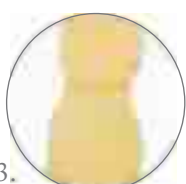

3.



4.

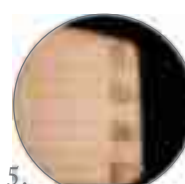

5.



6.

specifications

MAX. OPERATING WEIGHT	550 lbs.
WEIGHT	23-27 lbs
HEIGHT RANGE	21"-32"
WIDTH	28" x 21" x 28"
LENGTH	65"
TOP	Contour
PADDING (CFC FREE)	Semi-Firm, Plush, Firm
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Ultimate Access Ends
ACCESSORIES	QuickLock™ Face Rest Platform, Arm Hammock, Field Feet, Recess Fillers, and many more of our accessories

Bold - Standard • *Italic - Additional Charge*

- PATENTED INTEGRAHINGE™**
OAKWORKS® exclusive design integrates the attachment of the bracing system to offer incredible strength and lighter weight that no piano hinge can match.
- UNILOCK™**
Super quick, easy to use closure system. Won't rip your sheets.
- ERGONOMIC CONTOUR TOP**
Exceptional body mechanics for you and lighter weight for easy carrying. (Sample image shown in Saffron upholstery)
- COMPLEMENTARY SUSPENSION SYSTEM**
Designed to be the strongest system on the market while ensuring the best access to your client by allowing you to get your legs under the sides and the ends of the table, at any height!
- EXCEPTIONAL FRAME STRENGTH**
Frames are built on a computerized machine for accuracy and strength, utilizing box joints that are glued and screwed together.
- TELESCOPING LEGS**
Easy push button adjustment for quick height changes.

CS System™ (Complementary Suspension) : U.S. Patent #6,192,809
UniLock™: U.S. Patent # 7,472,656




*One great table...
One incredible value...*

ONE COMPLETE PACKAGE INCLUDES:

- 3" MULTI-LAYER PADDING
- ACCESS END & SHIATSU CABLES
- ADJUSTABLE FACE REST PLATFORM WITH AERO•CEL™ FACE REST CRESCENT
- ESSENTIAL CARRY CASE
- ARM HAMMOCK

specifications

MAX. OPERATING WEIGHT	500 lbs.
WIDTHS	30", 28"
LENGTH	73"
WEIGHT	25-32 lbs
HEIGHT RANGE	24" - 34"
TOP	Rectangular
PADDING (CFC FREE)	3" multi-layer padding
FABRIC (PVC FREE)	8 TerraTouch™ colors
STANDARD	Access End Panel, Shiatsu Cables
Bold - Standard • <i>Italic - Additional Charge</i>	

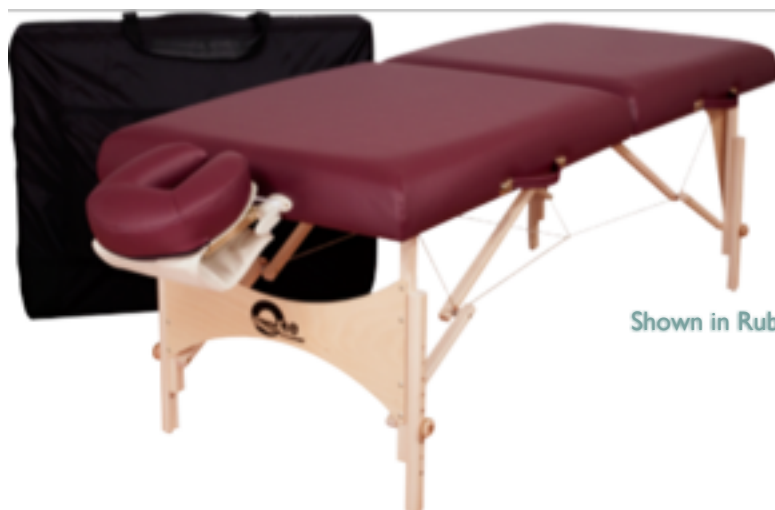

Shown in Ruby

Details Make the Difference



- 1. CABLE LOCK SYSTEM**
Our system uses stronger cables than anyone else, designed not to stretch or fail even under the heaviest loads.



- 2. SUPPORT STRUCTURE**
Engineered Leg extension system—designed to need only one knob instead of two and offer more strength than any other on the market.



- 3. UNDERSTRUCTURE**
This IntegraHinge™ design offers incredible strength and light weight that no piano hinge can match. Squeak-free and rigid.



- 4. HANDLES**
Reinforced, padded and T-nutted to the frame for added durability.



- 5. UNILOCK™ SYSTEM**
Super quick, easy to use system folds out of the way for shiatsu and won't rip your sheets.

UniLock™: U.S. Patent # 7,472,656

FABRIC OPTIONS



professional line MONEY SAVING

ESSENTIAL PACKAGE

INCLUDES:

- Essential Carry Case
- QuickLock™ face rest platform
- Aero•Cel™ face rest crescent
- 6" Semi-round Bolster

Package available for: Nova™, WellSpring™, Aurora™, Advanta™ Wood or Aluminum



PACKAGE NOT AVAILABLE WITH ADVANTA TABLE

Aurora™ Essential Package
shown in Bluegrass

TRAVELER PACKAGE

INCLUDES

- Professional Carry Case
- QuickLock™ face rest platform
- Aero•Cel™ face rest crescent
- 6" Semi-round Bolster
- Arm Hammock
- Table Cart

Package available for: Nova™, WellSpring™, Aurora™, Advanta™ Wood or Aluminum



WellSpring™ Traveler Package
shown in Heron

ACCESSORY PACKAGES

See pages 19-22 for accessory details.

PROFESSIONAL PACKAGE

INCLUDES

- Professional Carry Case
- QuickLock™ face rest platform
- Boiance™ face rest crescent
- 6" Semi-round Bolster
- Arm Hammock

Package available for: Nova™, WellSpring™, Aurora™, Advanta™ Wood or Aluminum



Nova™ Professional Package
shown in Espresso

ULTIMATE PACKAGE

INCLUDES

- Professional Carry Case
- QuickLock™ face rest platform
- Boiance™ face rest crescent
- 8" Fluffy Bolster
- Arm Hammock
- Flannel Sheet Set
- Premium Heated Table Fleece
- Premium Swivel Stool

Package available for: Nova™, WellSpring™, Aurora™, Advanta™ Wood or Aluminum



Alliance™ Wood Ultimate Package
shown in Sage

PORTAL PRO[®]

PACKAGE

professional line

The patented PortalPro[®] has been the undisputed favorite massage chair for over 20 years. This lightweight chair offers superior ergonomics, flexibility and strength. With a working weight of 350 lbs, the PortalPro[®] easily adjusts to fit any size client.

specifications

MAX. OPERATING WEIGHT	350 lbs.
WEIGHT	19 lbs
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Carry Case with strap, Instructional DVD
FRAME	Black Powder-Coated Aluminum Frame
ACCESSORIES	<i>Sternum Pad</i>

Bold - Standard • *Italic - Additional Charge*

QuickLock™ Face Rest: U.S. Patent #5,177,823
Portal Pro[®] : U.S. Patent #5,401,078



Package Includes:
Portal Pro[®] Chair,
Carry Case w/strap,
& Instructional DVD
(Shown in Sage)

- 1. STERNUM PAD OPTION**
Ultimate support for large breasted women, mastectomy patients, post-partum and pregnant clients.
- 2. ERGONOMIC DESIGN**
The chair instantly feels comfortable to people of all sizes. It features very high quality firm foam in the seat and knee pads for great comfort and longevity and utilizes softer foam in the chest and armrest for just the right support.
- 3. ACCESSIBILITY AND ADJUSTABILITY**
Ultimate access to your client's neck, shoulders, arms, and back. Simply adjust to support any size and any shape body. Match your working height and relieve strain on your client's lumbar spine.
- 4. WHEELS**
Out of the way when you are working, then instantly engaged for easy and fast traveling.

PORTAL LIGHT



Shown in Ruby with black trim

Weighing just over 15 lbs. the OAKWORKS® Portal Light will revolutionize the massage chair industry. Made with aircraft aluminum, no other light weight chair on the market offers the strength, ease of set up and client comfort of the new Portal Light.



convenient carry strap & carry case



personal item storage shelf

specifications

MAX. OPERATING WEIGHT	300 lbs.
WEIGHT	15 lbs
FABRIC (PVC FREE)	8 TerraTouch™ colors With or without black trim
PORTAL LIGHT PACKAGE INCLUDES	QuickLock™ Face Rest Platform & Aero•Cel™ Face Rest Crescent, Carry Strap & Case, Instructional DVD
ALUMINUM POWDER COATED FRAME	Black
ACCESSORIES	Sternum Pad
Bold - Standard <i>Italic - Additional Charge</i>	

QuickLock™ Face Rest: U.S. Patent #5,177,823

QUALITY FOAM & FABRIC

Oakworks professional TerraTouch™ 100% polyurethane fabric with a 5 year guarantee and high quality US foam.

ERGONOMIC DESIGN

The Portal Light is ergonomically engineered to be comfortable for virtually all clients.

EASE OF USE

Easy set up with very little adjusting means you can use your time working on clients.

LIGHT WEIGHT & DURABLE

Made of aircraft quality aluminum, the Portal Light is very durable, safe and stable. At only 15 lbs, you will want to take it everywhere!

FABRIC OPTIONS



DESKTOP PORTAL™



Shown in Sage

specifications

WEIGHT	11 lbs.
FABRIC (PVC FREE)	TerraTouch™
FRAME	Black Powder-Coated Aluminum Frame
STANDARD	Carry Case with strap
ACCESSORIES	Sternum Pad
Bold - Standard • <i>Italic - Additional Charge</i>	

Expand your practice to include everyone, even clients who are bed-ridden, desk-bound or just really busy. Corporate offices, hospitals, nursing homes, hotels... all within your reach. Thanks to the Desktop Portal™'s light weight and superior ergonomics, the opportunities are endless...



STORABLE MAT

FEATURES:

- INCREDIBLY STABLE DESIGN
- 33", 35", EVEN 40" WIDTHS
- 16"-25" HEIGHT RANGE



Shown in Sapphire with optional stabilizing feet

CS System™ (Complementary Suspension):
U.S. Patent #6,192,809

UniLock™: U.S. Patent # 7,472,656

The most stable portable table on the market! Perfect for all forms of structural work, the Storable Mat allows you to use good body mechanics while working on a surface as stable as the floor with its 16" to 25" height range. Oakworks patented Complementary Suspension System is the key to the incredible stiffness. Optional Stabilizing Feet permit super stable work right on the edge of the table.

specifications

MAX. OPERATING WEIGHT	550 lbs.
TABLE WEIGHT	40-55 lbs (depending on width)
WIDTHS	33", 35", 40"
LENGTH	78"
HEIGHT RANGE	16" - 25"
TOP	Rectangular
PADDING (CFC FREE)	2.5" Semi-Firm <i>1.5" Firm Response</i>
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
OPTIONS	Stabilizing Feet (4 feet included with 40" width table)
Bold - Standard • <i>Italic - Additional Charge</i> (except padding & 40" width)	



STABILIZING FEET

For maximum stability each 40" Storable Mat comes standard with stabilizing feet, which easily fold out of the way for storage.



PROLUXE™

professional line

The ProLux™ electric lift table provides exceptional value with its Dual X Brace design providing unmatched stability and wide range of tops.



Foot control included!

Shown with Flat top & Ocean fabric and optional QuickLock™ face rest platform & Boiance™ face rest crescent.



Shown with Lift-Assist Backrest top & Orchid fabric and optional Large Storage Shelf, Hot Towel Cabi Shelf and Power Outlet


TOP & PADDING OPTIONS | 29" or 31" widths | 73" length



Electric Salon
4" Aero•Cel™ standard,
6" Aero•Cel™



Lift-Assist Salon
4" Aero•Cel™ standard,
6" Aero•Cel™

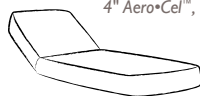

Lift-Assist Backrest
3" Plush standard,
4" Aero•Cel™, 6" Aero•Cel™



Flat
3" Plush standard,
4" Aero•Cel™, 6" Aero•Cel™



ProLux™ shown here with accessories from Comphy.

specifications

MAX. OPERATING WEIGHT	550 lbs.
HEIGHT RANGE (with 4" Aero•Cel™ padding)	18"-35" (46-89cm) flat top, 21"-38" (53-97cm) backrest & salon tops
HEIGHT RANGE (with optional shelves & 4" Aero•Cel™ padding)	25"-36" (64-91cm) flat top 27"-39" (69-99cm) backrest top 27"-39" (69-99cm) salon top
WIDTHS	29", 31"
LENGTH	73"
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	Multi-function foot control
OPTIONS	Twin locking casters, Hand Control, Large Shelf, Cabi Shelf, Power Strip
ACCESSORIES	QuickLock™ face rest platform
Bold - Standard • <i>Italic - Additional Charge</i>	



Large Storage Shelf (all tops)
8.75" H x 19.75" W x 13" D



Hot Towel Cabi Shelf (varies by top)
Flat Top - Fits up to 13" H x 19.75" W x 9" D
Backrest & Salon Tops - Fits up to
12" H x 19.75" W x 9" D



Power Strip
Included with every Hot Towel Cabi Shelf
allows the user to plug everything in here,
and have only one cord to plug in to your
power source.



PROLUXE™ CONVERTIBLE

converts in seconds! ...then back to a portable!

ELECTRIC HEIGHT ADJUSTMENT

Quickly change the height of your table during treatments without disturbing your client.

USE YOUR PORTABLE TABLE


Use the table you already have and convert it into a lift table! Works with portable tables that can be laid flat and are 27"-33" wide, 72"-75" long.

EASE OF USE

Converts in less than a minute! Use the lift base in your office and take your portable for outcalls.

CONTROLS INCLUDED

Choice of hand or foot control.

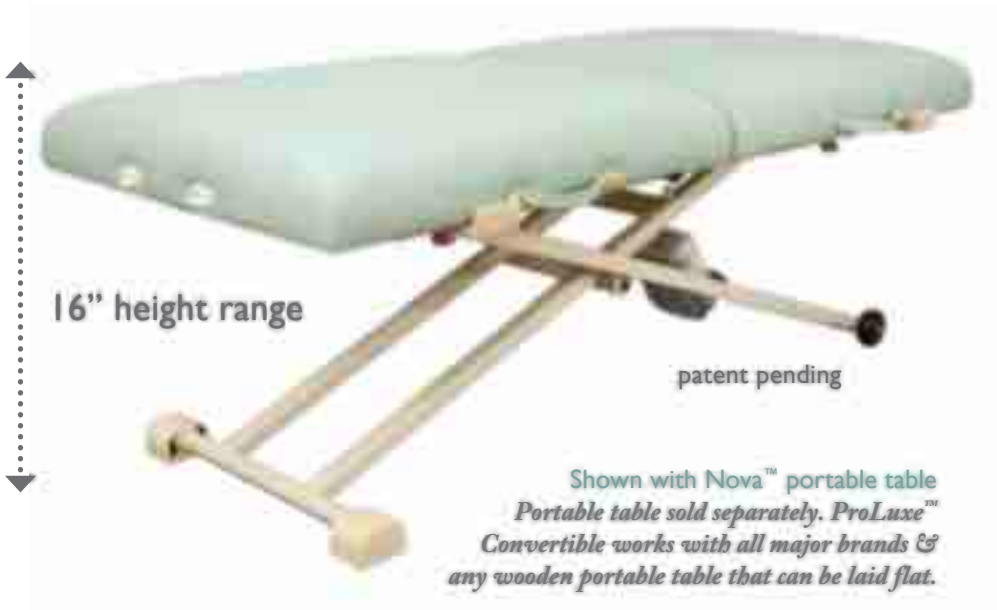

specifications



MAX. OPERATING WEIGHT	Lifting capacity is 550 lbs. Weight is defined by the table that is mounted to the unit.
DIMENSIONS	29" W x 66" L Fits any portable table that is 27"-33" wide, 72"-75" long
PRODUCT WEIGHT	115 lbs.
CONTROL OPTIONS	Hand Control or Foot Control
OPTIONS (additional charge)	Extra controls, Locking Twin Casters




.....> your portable table> converted to a stationary lift table> and back again>



RECOMMENDED TABLES:
Nova™, Alliance™ Wood, Aurora™, One
or any brand massage table that can be laid flat!

you deserve
the **PROLUXE™**
convertible
patent pending

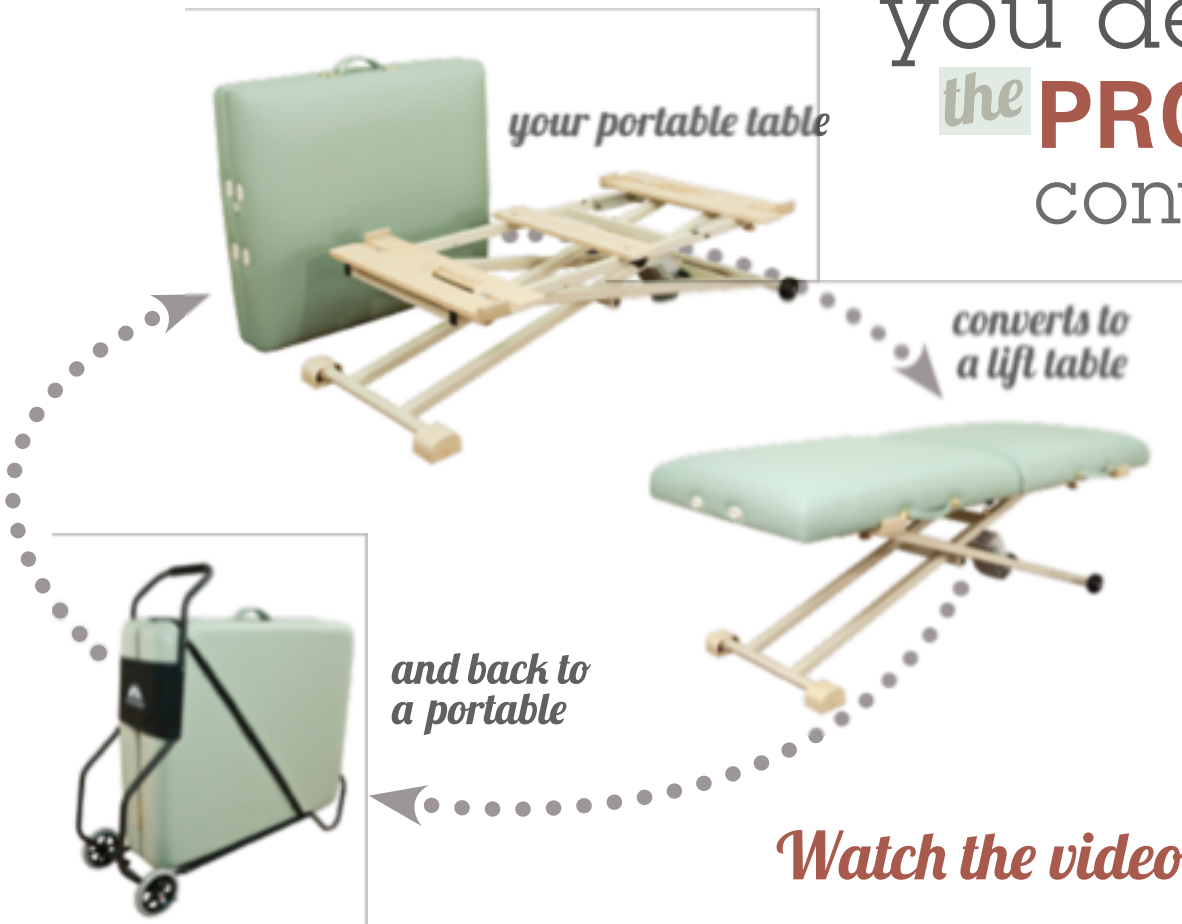

Watch the video!

*tested
& approved!*

OAKWORKS®
EDUCATIONAL SERIES

Working together to promote massage education



“I have long felt that an electric lift table is not a luxury, but an essential tool for your practice. Your clients will benefit from improved care and your business will benefit as a result!”

Joe Muscolino
www.LearnMuscles.com



“The Oakworks ProLuxe™ brings together the best of both worlds. Client comfort and therapist functionality. Customize the height of your table for every client, every session!”

Susan Salvo
www.LAMassageSchool.com



“Finally, a massage table lift that is affordable and easy to add to your practice. The lift allows me the freedom to move from my table shiatsu practice to pregnancy massage to foot reflexology sessions. Versatility at its best!”

Claire Marie Miller
www.ClaireMarieMiller.com



“You deserve a hi-lo table as a therapist. You need to save your own body if you want to continue working. Your table will save your body!”

Taya Countryman
www.StructuralReliefTherapy.com

PT400M *professional line*

The PT400M is the most advanced massage table on the market. It provides excellent access to your clients while providing comfortable support. The powered Center break in prone opens up their back muscles and relieves pressure in the lower back. In supine it takes the pressure off their lower back and creates a gentle abdominal stretch. The powered Center Break also allows you to continuously stretch muscles in side lying, prone, and supine by slowly raising the center as the session progresses.

- **Advanced design for the best ergonomics**
- **Powered Center Break and Height**
- **Qualifies for ADA Tax Credit**



The PT400M comes with QuickLock™ Face Rest Platform, Boiance™ Face Cradle Pillow and Two-function Foot Control.



PRONE

The PT400M relaxes and stretches the lower back naturally. Shown with Extension Pad.



SUPINE

The PT400M supports the lower back and gently stretches the abdominals.



SIDE ARM RESTS

Side arm rests allow you to support your client and take advantage of the narrow 22" width. They pivot in and out and can be easily removed.



POWERED CENTER BREAK



Watch the video!



EXTENSION PAD


The extension pad can be used at the head or foot end, and as an arm shelf with the arm shelf bracket accessory.



specifications

MAX. OPERATING WEIGHT	550 lbs.
HEIGHT RANGE	17"-34" (44-87cm)
WIDTH	29.5"
LENGTH	77" (including head rest)
PADDING	2.5" Comfort Foam™ <i>1.75" Firm Response, 4" Aero•Cel™</i>
FABRIC (PVC FREE)	TerraTouch™ choose from over 20 colors!
STANDARD	QuickLock™ Face Rest Platform & Boiance™ crescent & Foot Control
OPTIONS	<i>Hand Control, Second Foot Control, Locking Twin Casters</i>
ACCESSORIES	<i>PT400M Rotating Arm Rest (set of 2), PT400M Extension Pad, PT400M Shelf Bracket,</i>

Bold - Standard • *Italic - Additional Charge*

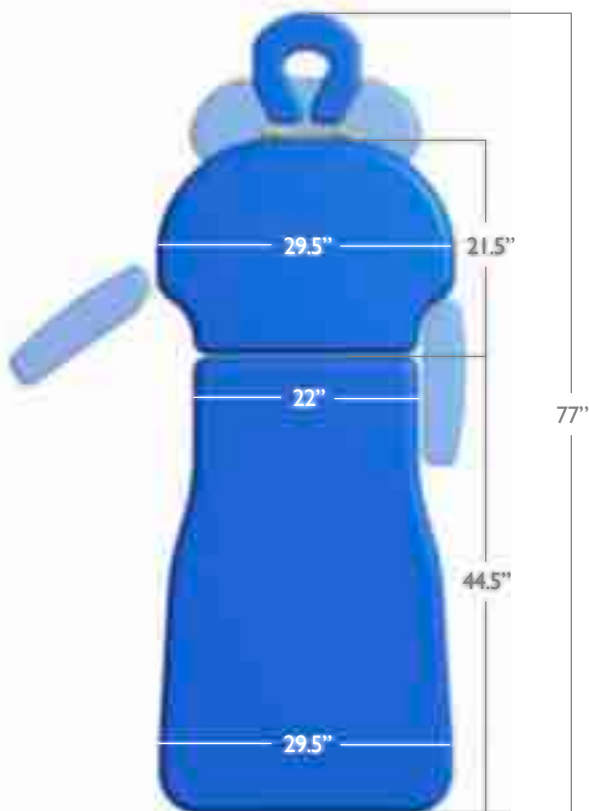

CLINICIAN™ *professional line*

specifications

MAX. OPERATING WEIGHTS	**550 lbs., *500 lbs.
HEIGHT RANGE	24"-36"
WIDTHS	29", 31"
LENGTH	73"
FABRIC (PVC FREE)	TerraTouch™
STANDARD	Storage shelf
OPTIONS	Wood finishes, Optional cabinet w/ door designs and our full line of accessories
ACCESSORIES	QuickLock™ face rest platform
Bold - Standard • Italic - Additional Charge	



OAKWORKS® Clinician™ table line provides incredible versatility and comfort. Unlike many spa treatment tables, OAKWORKS®' proprietary Clinician™ design incorporates a higher quality and unique construction details that make it the most durable, reliable table. Loaded with features able to withstand the test of time, 3 top options and 3 lift systems, it's the simple choice for your stationary table!

TOP & PADDING OPTIONS | 29" or 31" widths | 73" length



Lift-Assist Salon
4" Aero+Cel™ standard,
6" Aero+Cel™



Lift-Assist Backrest
3" Plush standard,
4" Aero+Cel™, 6" Aero+Cel™



Flat
3" Plush standard,
4" Aero+Cel™, 6" Aero+Cel™

MAX OPERATING WEIGHTS

*Clinician™ with Adjustable Lift supports 550 lbs.

**Clinician™ with Manual/Hydraulic or Electric/Hydraulic Lift supports 500 lbs.

ADJUSTABLE LIFT

Wooden engraved knobs for easy height adjustments.



Manual leg height adjustment



550 lb. Weight Capacity

Shown with Flat top, Heron fabric, Natural finish & storage shelf and optional QuickLock™ face rest platform & Boiance™ face rest crescent.

WOOD FINISHES

Available on all Clinician™ tables! Ask for details on our stains and lacquers to complete your table's look.

STAINS




Cherry



Golden Oak



Walnut

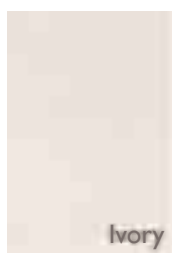

Ivory

LACQUERS




Cabernet



Ebony



Chocolate

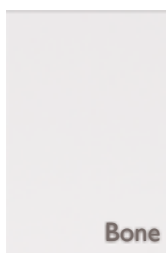

Bone

MANUAL/HYDRAULIC LIFT

This hybrid system replaces the electric motor of the electric/hydraulic system with a fold away crank handle that easily lifts or lowers your client to the desired height. A great choice where electricity is problematic.



Manual crank height adjustment



Folds out of the way

500 lb. Weight Capacity

Shown with Lift-Assist backrest top & Seafoam fabric and optional Cabinet with Standard One Piece Frosted Acrylic doors & Walnut stain.

ELECTRIC/HYDRAULIC LIFT

Quick, quiet adjustment without disturbing your client.

Foot control is included.

LIFT-ASSIST SALON & BACKREST TOPS

The powerful gas cylinders are controlled by easy to reach finger controls.



Foot control included!

500 lb. Weight Capacity

Shown with Lift-Assist salon top, Taupe fabric & Natural finish and optional Cabinet with Standard Frame and Panel Wood doors, Universal Table Extender & Adjustable Side Arm Rests.



CABINET OPTIONS

Available on all Clinician tables! Ask for details on all of our Standard & Signature door designs & finishes.



Standard Frame and Panel Wood doors in Natural finish



Signature Sierra Mist One Piece Solid Surface doors in Cherry Stain



Boiance™

FACE REST CRESCENT

Watch the video!




The only face rest crescent with built in water spheres for cushion and support!

the
MOST COMFORTABLE
FACE REST CRESCENT
ever

PATENT PENDING WATER SPHERE TECHNOLOGY

introducing

Boiance™ float

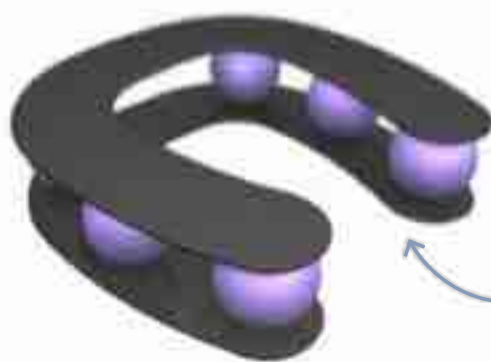

shown with the Boiance™
face rest crescent

- Incredible freedom of movement in all directions
- Eliminate compression of the cervical spine
- Works with most major brands of face cradles & pads
- Patent Pending Water Sphere Technology
- Just insert between cradle and pad!



*Watch
the video!*

The Boiance Float device is excellently designed. It allows the head and neck to move naturally as the massage proceeds. My clients say that it is the most comfortable their neck has ever been on a massage table!

Dr. Ben Benjamin

The Boiance Float comes with a 30 Day No Risk, 100% money back guarantee! If you're not satisfied with this product, simply return it within 30 days from the date of purchase for a full refund of the purchase price and shipping charges. *Available in Black only.*



ACCESSORIES

Visit www.massagetables.com to see
our entire line of bolsters & accessories!

OAKWORKS® BOLSTERS

100% polyurethane fabric, no PVC, pull handles at both ends.



6" Half-round
3"H x 26"L

8" Half-round
4"H x 26"L

6" Semi-round
4.5"H x 26"L

Most Popular!

8" Semi-round
6.5"H x 26"L

6" Round
6"D x 26"L

8" Round
8"D x 26"L

AERO•CEL™ FACE REST CRESCENT

Soft and supportive, this versatile face rest pad adjusts to support any size or shape face for prone comfort. Turn it around and create a great cervical pillow for supine treatments. It is the standard pad included with every face rest.




QUICKLOCK™ FACE REST PLATFORM

Our best-selling face rest platform is an invaluable support option, providing superior access with optimum comfort. The double articulating adjustability ensures a perfect fit for all body types and even helps with positioning clients with cervical injuries. (Available in Black or Vanilla)



ADJUSTABLE FACE REST PLATFORM

This face rest adjusts in height and angle using a single cam lock.

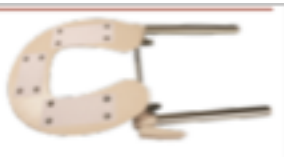

BUTTERFLY CUSHION

The adjustability afforded by the air system ensures proper positioning and comfort. The two chamber design, with asymmetrical support regions, offers a wide range of positioning options.



Watch the video!

FLUFFY BOLSTER

OAKWORKS® 8" Fluffy Bolster is the perfect positioning tool offering cradling support for the hips, knees and ankles. Filled with an ultra-soft material and covered with earth friendly PVC free TerraTouch™ fabric. 26"L with 8" Diameter



ADJUSTABLE 8" AIR BOLSTER

The patented 8" round air bolster features twin air chambers that allow you to customize the size of the bolster to conform perfectly to your client and offer incredible comfort compared to foam or poly filled bolsters. 26"L with 8" Diameter



BODY CURVE™

Prevents the client's neck from being hyper-extended during comprehensive treatments like facials, waxing and massage services. The angled sides offer accessibility without sacrificing support. The Body Curve™ can also be used for lower back support, or to support the wrist or shoulder. 11"L x 5"W x 3" H



WEDGE BOLSTER

Use a wedge bolster along with our QuickLock™ face rest platform to create a quick and easy backrest for any portable or stationary massage table. 24"W x 10"H x 24"L with a 36° Incline.



FITTED DISPOSABLE

FACE REST COVER

50 per box. Provides a fitted disposable sanitary covering at an irresistible price.



FLAT DISPOSABLE

FACE REST COVER

100 per package. Super soft and disposable, this flat cover provides a sanitary covering for just pennies per treatment.



KING SIZE PILLOW COVER

18" W x 34" L. Fits king size pillow 20" x 36" or two standard pillows. Available in all TerraTouch™ colors. Includes pull handle at one end and zipper closure.



FITTED TABLE PROTECTOR

Now you can easily protect your massage table with this PVC free TerraTouch™ fitted table protector. Fits tables up to 31" wide.




ESSENTIAL CARRY CASE

A sturdy, basic case with a large zippered pocket and multiple handles. It gives you the table protection and carry convenience you need at an affordable price. Sizes: L, XL



PROFESSIONAL CARRY CASE

A large 19"x 24" storage pocket, padded bottom for added protection, a non-slip shoulder strap, convenient handles and dual-slide zippers take you comfortably down the road, while the durable case protects your table in style. Sizes: M, L, XL, XXL & XXXL

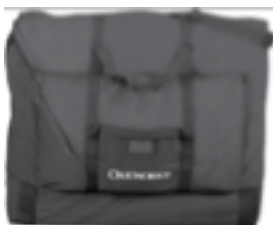

WOODEN PORTABLE ARM REST

Available on the Nova™ and Aurora™. Hinged to the table. Folds out of the way when not in use. Not available on table with height range of 18" - 26". Face Rest required.




ARM HAMMOCK

Offer your clients optimal comfort without spending a fortune. The OAKWORKS® Arm Hammock is the "must have" massage table accessory. This unique product will comfortably and easily support your client's arms when in the prone position. Available in Black or Vanilla. Face rest required.




POLY-COTTON SHEETS

Includes: Fitted Sheet, Flat Sheet & Face Rest Cover
Poly/cotton sheets wrinkle less, shrink less and last a very long time. Fits tables up to 31" wide. 200 Thread Count. Available in Chocolate.

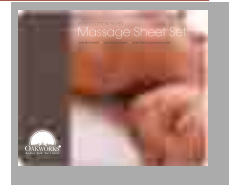

100% COTTON FLANNEL SHEETS

Includes: Fitted Sheet, Flat Sheet & Face Rest Cover
Machine washable, luxurious, unbleached, 100% cotton flannel sheets. Fits tables up to 31" wide. Available in Ivory.

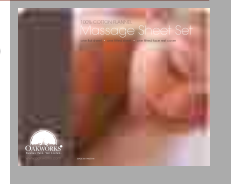

ESSENTIAL TABLE WARMER

The Essential Table Warmer adds a layer of warmth to your table, comforting your client and preparing them for a relaxing treatment. 30"W x 72" L



PREMIUM TABLE WARMER

Luxurious, fleece table warmer fits most tables and distributes therapeutic heat evenly while enveloping your clients in plush comfort. Programmable up to 99 minutes. 30" W x 72" L



TABLE & FACE REST FLEECE

Protects your table & face rest from exposure to oil, it's washable, and provides your clients with warmth and comfort. 28"W x 72" L



SIDE LYING POSITIONING SYSTEM



CAROLE OSBURNE
Author of Pre- and Perinatal Massage Therapy
& 2008 AMTA National Teacher of the Year

"I have been teaching pre- and perinatal massage for thirty years, and this system does it all! It is the first system designed specifically to comfortably support and stabilize every part of the body for side lying perfection."



GEORGE KOUSALEOS, BA, LMT, NCTMB
Founder and president of the CORE Institute School of Massage Therapy and Structural Bodywork in Tallahassee, Florida

"The SLP system is the finest quality bolstering system that I have ever used. It adds another level of professionalism to my Structural Integration and Myofascial Therapy practice."



BENNY VAUGHN LMT, ATC, NCTMB, CSCS
Licensed Massage Therapist since 1974. He is a Certified Athletic Trainer and a Certified Strength and Conditioning Specialist

"Placing clients and patients, in the most comfortable and least stressed position, when administering hands on, therapeutic intervention will magnify the best outcome. The ability to work with people in a side lying position, comfortably, will allow the therapist to explore more therapy options."



- ✓ Helps relax and stabilize your client to make your work easier
- ✓ Stretch soft tissues and mobilize joints with more ease
- ✓ Perfect for pregnancy massage, medical massage, sports massage
- ✓ 30 day money back guarantee
- ✓ Comes with free instructional DVD

No other system on the market is designed specifically to create ergonomically correct support for side lying positions. In development for almost a year with the help of top pregnancy instructors who are very focused on correctly positioning people in the most comfortable way, every detail was considered. Each pad has unique features that insure that the client is properly supported.



Watch the video!



STOOLS

These durable, comfortable stools support up to 300 lbs. and offer simple pneumatic height adjustments.



BASIC STOOL

Colors: Black, Heron or Ruby ONLY
Height Range: 18" - 24"



PREMIUM STOOL

Colors: All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 20" - 28"



PROFESSIONAL STOOL

All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 21" - 28"



PREMIUM STOOL with SADDLE SEAT

Colors: All TerraTouch™ Colors
Low Height: 19" - 25"
High Height: 22" - 29"



PROFESSIONAL STOOL with SADDLE SEAT

All TerraTouch™ Colors
Low Height: 19" - 25"
High Height: 22" - 29"



PREMIUM STOOL with BACKREST

Colors: All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 20" - 28"



PROFESSIONAL STOOL with BACKREST

All TerraTouch™ Colors
Height Range: 18" - 24"

TABLE CART

We designed the Table Cart specifically for Massage tables. Most on the market are luggage carts being sold as massage table carts. This design has soft support panels that cradle your table to insure the foam and fabric are not dented or ripped. We also balanced it perfectly, so when you are pulling your table around, it isn't pulling you down! Fits tables up to 35" wide. Available in black only.



TERRATOUGH™ FABRIC

TerraTouch™ Oakworks' PVC- free premium polyurethane fabric has excellent resistance to oil, abrasion, and tears. TerraTouch™ passes testing for Cytotoxicity, Skin Irritation and Sensitization as well as CA TB117 Fire Retardance testing. For disinfecting purposes you may use Protex, MadaCide, MadaWipe, Accell TB, Virox® (Canada), PDI Duper Sani-Cloth Germicidal wipe or a diluted 1/10 bleach solution.



PADDING

See product specifications for available padding on specific tables.

AERO•CEL™ PADDING



New technology in padding options now gives you 4" (6" available on stationary tables only) of more resilient support with a softness that conforms to the body's curves. Offering the most relaxing and nurturing environment possible, it's an investment that will return long years of "ahhhhs". Weighs 4-5 pounds more than semi-firm.



SEMI-FIRM PADDING



2 ¼" thick, 2 layer multi-density system is a great balance of lightweight and comfortable support. All portable table weights are based on this system.

FIRM PADDING



1 ½" thick, 2 layer multi-density system is our firmest option.

FIRM RESPONSE PADDING



1 ¾" thick, 2 layer multi-density system is designed specifically for physical therapy. Available only on the PT400M.

PLUSH PADDING



3" thick, 3 layer multi-density system creates the perfect balance of comfort, longevity, and light weight for portable tables.

COMFORT FOAM™ PADDING



2 ½" thick, 2 layer multi-density system is similar to Plush padding in feel, but is designed for the medical market and higher expected usage rates. Available only on the PT400M.



OAKWORKS® EDUCATIONAL SERIES

Working together to promote massage education



Free

ONLINE VIDEO LIBRARY
FROM TOP PROFESSIONALS

- Demonstrations using Oakworks® products
- Learn techniques from educators in the field
- Take your career to the next level!

Scan and watch now

YouTube

www.youtube.com/oakworksinc



Educators

*over 200 videos!
and growing...*

Carole Osborne, CMT
Joseph E. Muscolino, DC
Walt Fritz, PT
Susan Salvo, LMT
Judy Koch, CMT
Allison Ishman, LMT, SMT, PFT
Debbie Roberts, LMT
Taya Countryman, LMP
Paul Lewis, RMT
Jeffrey Haggquist
Roger McNear
Claire Marie Miller, LMBT
Samuel Blechman, LMT
Annie Powell

Ben Benjamin, Ph.D
Aaron Mattes, LMT
Ralph Stephens, LMT
Kathleen Gramzay, LMT
Irene Diamond, RT
Kate Jordan, NCBTMB
Dr. Bob Duvall
Nancy Porambo, MS, CNMT, NCTMB
Gloria Coppola, LMBT
Patrick Ingrassia, LMT
Katie Adams, LMT, CNMT, CMTPT
Stewart Wild, LMT, CNMT, CMTPT
Pat Mayrhofer, CMT
and many more...



717.235.6807
www.massagetables.com