

e-xam



Technical data:

e-xam LED 2.5 V

5.000-5.500 Kelvin

180.000 Lux max.

field of light: approx. 11 cm diameter
(distance 15 cm)

e-xam Xenon 2.5 V

3.500-4.000 Kelvin

60.000 Lux max.

field of light: approx. 14 cm diameter
(distance 15 cm)

Batteries: 2 x Alkaline AA, 1,5 V



e-xam

Diagnostic pen-light

Photobiologically validated high-performance pen-lights for medical diagnosis of throat and pupil.

Available with LED or Xenon illumination.

- Photobiologically validated for lighting of pupils according to EN 62471:2008
- Pre-focused and even field of light.
- Bright and natural white light by high-performance LED.
- LED burns at least 20.000 hours.
- Low maintenance costs thanks to long LED and battery service lives.
- ABS casing with metal clip, disinfectable.
- easy to use with on/off switch.
- including blade holder, batteries and nylon case.
- Designed for durability, comes with a 2-year warranty including LED.

Ordering information:


e-xam

No. 5130-01 black, XL 2.5 V, with tongue blade holder

No. 5130-02 white, XL 2.5 V, with tongue blade holder

No. 5131-01 black, LED 2.5 V, with tongue blade holder

No. 5131-02 white, LED 2.5 V, with tongue blade holder

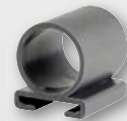

2,5 V LED for e-xam

No. 12320, pack of 1

XL 2.5 V Xenon bulb for e-xam

No. 11178, pack of 6

fortelux N



Technical data:

3 V LED

5,500 Kelvin

24,000 Lux max.

Field of view approx. 50 mm diameter
(distance 15 cm)

Batteries: 2 x Alkaline AAA, 1,5 V

fortelux N

Diagnostic pen-light

Economy pen-light for medical diagnosis with LED illumination

- LED illumination close to natural daylight.
- Pre-focused field of light.
- Affordable price.
- Aluminium outer casing with metal clip, disinfectable.
- Simple on/off function through contact with the metal clip.
- Including batteries.
- Optional blade holder, laryngeal mirror and blue filter.
- Comes with a 2-year warranty.



Ordering information:

fortelux N, pack of 1

No. 5075-526 black

No. 5074-526 silver

No. 5077-526 red

No. 5071-526 blue

No. 5078-526 white

No. 5072-526 green

fortelux N, pack of 6

No. 5070 assorted colours

No. 5075 black

No. 5074 silver

No. 5077 red

No. 5071 blue

No. 5078 white

No. 5072 green

ri-light set (fortelux N, blade holder, laryngeal mirror, blue filter)

No. 5090 white



3 V LED for fortelux N

No. 12322, pack of 1